

The Spectroscopy Specialists

Quartz and Glass Cells for:





about Starna Cells, Inc....

The origin of the optical skills available from **Starna**, formed in 1971, can be traced back to the earlier part of the 20th century. Utilizing their optical expertise, the founders of the original company manufactured specialized optics and, during the early 1950s, had already developed the technique necessary for manufacturing the first *fully fused* form of spectrophotometer cells.

As a wholly owned subsidiary of the **Starna** international group of companies, Starna Cells provides many instrument manufacturers with specialized optical components and other private labeled products, including the Starna[®] brand of high quality spectrophotometer cells, accessories and the extensive range of NIST Traceable Reference Materials.

Starna products are distributed worldwide through the Starna network to instrument manufacturers, distributors and direct to end users. A dedicated technical applications and development team is ready to answer questions.

Cell Construction Specifications

All the cells Starna manufactures, unless specifically designed otherwise, are assembled using a fully fused method. This technique, pioneered and perfected by Starna, ensures that cells are fused into a single homogeneous unit by heat alone, using no intermediate adhesives. The cells achieve maximum physical strength as well as being unaffected by solvents.

All cells are carefully annealed to remove any residual strain left from the fusing process and, with few exceptions, can be used safely with pressure differentials of up to three atmospheres.

General Specifications:

Windows parallel to: Window flatness to: Window polish: Window thickness:	better than 3 minutes of arc less than 4 Newton fringes 60/40 scratch/dig 1.25mm				
Material-	Path lengths	Tolerance			
Glass	up to 20mm	+/- 0.1mm			
Glass	30 to 100mm	+/- 0.2mm			
Special Optical Glass	up to 20mm	+/- 0.01mm			
Special Optical Glass	30 to 100mm	+/- 0.02mm			
Quartz	up to 0.05mm	+/- 0.002mm			
Quartz	0.1 to 0.4mm	+/- 0.005mm			
Quartz	0.5 to 20mm	+/- 0.01mm			
Quartz	40 to 100mm	+/- 0.02mm			

The cells can be used with most solvents and acidic solutions. Acids such as HF should be avoided as they will attack the quartz itself while most other acids will not affect the cells. Strongly basic solutions (pH 9.0 and above) will etch the surface of the windows and shorten the useful life of the cells.

Note: Dimensions and specifications may change without notice.

Material Specifications

Starna offers five window materials, Optical Glass, Special Optical Glass and Pyrex[®] for the visible range, Spectrosil[®] Quartz for the far UV range and Infrasil[®] Quartz or equivalent for the near infrared. If a material required is not shown in this catalog, please contact us for availability. All materials used to construct the cells are suitable for use in the following wavelength range:

Optical Glass	334 through 2500 nm
Special Optical Glass	320 through 2500 nm
Pyrex [®]	320 through 2500 nm
Spectrosil [®] Quartz	170 through 2700 nm
Infrasil [®]	220 through 3800 nm

For fluorescent applications, Spectrosil[®] is recommended as it does not exhibit any background fluorescence. Some other materials, especially glass and lower grades of quartz may have some background fluorescence. Standard window thickness is 1.25mm, polished to a flatness tolerance of better than 4 Newton Fringes per centimeter in the viewing area. They are typically flat to better than 1 micron (0.001mm) over the whole window area. The scratch and dig specification for surface polish is 60/40. The meticulous care taken in the preparation and construction of regular quartz fluorescent cells allows for normal tolerances to be sufficiently stringent for use in laser applications.

Cell Matching

The high degree of accuracy maintained during production ensures a standard path length tolerance of +/- 0.01mm. The tolerance maintained for parallelism of the windows is better than 3 minutes of arc, therefore quartz cells vary little in transmission values. Slight differences in the transmission of new glass cells are due to variations in raw material transmission characteristics. The transmission of matched cells is measured and each cell is given a match code. These codes are only of real value when comparing new cells because the transmission characteristics may change as surface contamination or deterioration occurs during normal use. Thus, a new cell of a particular match code will not necessarily match an older cell of the same match code.


Cells which are commonly used (eg. 1-Q-10) are normally supplied from stock in unlimited quantities of the same match code. The less commonly used cells can usually be supplied in matched sets of two or four. In the case of Tablet Dissolution flow cells, matched sets of eight or more can be supplied.

Window	Matching	Measured at
Material	Tolerance	Wavelength
Optical Glass	0.5 %	350 nm
Special Optical Glass	1.0 %	320 nm
Spectrosil Quartz	1.5 %	200 nm
Infrasil Quartz	1.5 %	240 nm

'Z' Dimension for small volumes

The 'Z' dimension is the distance from the base of a cell to the center of the sample chamber window. The 'Z' dimension is very important for small volume cells of any design, where the sample compartment cross section dimension is very small. The correct 'Z' dimension should be added to the part number for small volume cells

'Z' Dimension per Instrument								
Manufacturer:	'Z' Dimension:							
Agilent [®]	15 mm							
Beckman [®]	8.5 mm							
Bio-Rad [®]	8.5 mm							
Eppendorf [®]	8.5 mm							
GBC ^(R)	15 mm							
Hewlett Packard [®]	15 mm							
	ries by instrument							
Jasco [®]	11 mm							
Ocean Optics [®]	15 mm							
Perkin-Elmer [®]	15 mm							
Pharmacia [®]	15 mm							
Shimadzu®	15 mm							
StellarNet [®]	15 mm							
Thermo Spectronic [®]	8.5 and 15 mm							
Turner [®]	8.5 mm							
Varian [®]	20 mm							

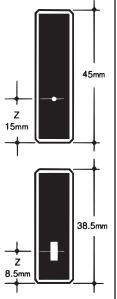

How to determine 'Z' Dimension:

1) Cut a piece of paper 12mm x 100mm

2) Punch a hole with a pen 15mm from one end and 8.5mm from the other end

3) Put your instrument in %T at about 535nm.

4) Place the paper in the cell holder and see which hole transmits light through the hole.

5) If neither transmits light, cut another piece of paper and try other dimensions from the end of the paper.



Contents:

Spectrophotometer Cells

Speciro	photometer Cells	
	Closed screw cap	12, 13
	Cylindrical	22 - 25
	Demountable, short path	14
	Demountable, flow	15
	Dual pathlength	30
	Flow cells	15 - 21
	Graded seal attached	35
	Low head space	11
	Micro volume	8, 9
	Mixing Cells	31, 37
	Semi-micro volume	6, 7
	Septum cap cells	12, 13
	Standard rectangular	, 4, 5
	Stirring cells	37
	Sub-micro volume	10
	Tablet dissolution	17 - 19
	Tube attached to cells	35
	Water jacket cells	38
Fluoron	neter Cells	_
	Closed screw cap	31
	Dual pathlength	30
	Flow cells	32, 33
	Graded seal attached	35
	Micro volume	26, 27
	Micro, square	28
	Semi-micro volume	26, 27
	Septum cap cells	31
	Standard rectangular	26, 27
	Stirring cells	37
	Sub-micro volume	29
	Tandem or Divided	31
	Triangular	30
	Tube attached to cells	35
Special	Cells and Accessor	ride
opeciai		
	Cell adaptors 7, 14, 15, 24 Colorimeter cells	
		36
	Laser cells	34
	Magnetic Stirrer	37
	Microscope cells	16
Validati	on sets	
	Fluoresence reference set	33

Standard Rectangular

1mm


2mm

5mm

Standard rectangular cells or Macro cells, Type 1, are the most common cells used in analytical chemistry. The exterior dimensions allow the cells to be used with virtually every spectrophotometer and photometer instrument currently in use.

These cells have an exterior height of 45mm, two clear windows and an open top. The walls of the cell and the bottom are smooth on the inside. On the outside, the light path windows are clear while the adjacent windows and the bottom are frosted or grey. Typically, windows are 1.25mm thick and each cell in this range is supplied with a PTFE cover.

10mm

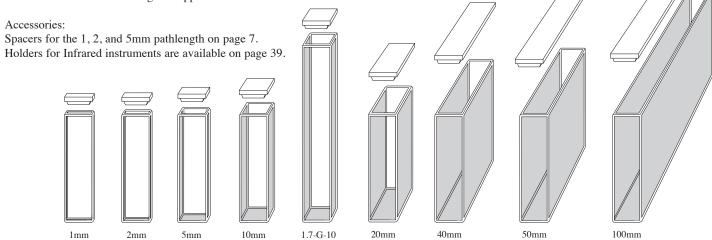

20mm

1.7-G-10

50mm

40mm

Catalog Number	Path Length mm		xterior, m Length		Interi Width	or, mm Length	Nominal Vol. ml	Polished Windows	Price per cell
Optical Glass	s windows, Useabl	e range: 3	334 to 250	0 nm					
1-G-1	1	12.5	3.5	45	10	1	0.400	2	\$64.00
1-G-2	2	12.5	4.5	45	10	2	0.700	2	64.00
1-G-5	5	12.5	7.5	45	10	5	1.700	2	58.00
1-G-10	10	12.5	12.5	45	10	10	3.500	2	42.00
1.7-G-10	10	12.5	12.5	70	10	10	6.500	2	90.00
1-G-20	20	12.5	22.5	45	10	20	7.000	2	64.00
1-G-40	40	12.5	42.5	45	10	40	14.000	2	81.00
1-G-50	50	12.5	52.5	45	9.5	50	17.500	2	111.00
1-G-100	100	12.5	102.5	45	9	100	35.000	2	189.00
Special Optic	al Glass windows,	Useable	range: 32	0 to 2500 n	ım				
1-SOG-10	10	12.5	12.5	45	10	10	3.500	2	\$58.00
Spectrosil [®] I	Far UV Quartz wir	ndows, Us	seable ran	ge: 170 to	2700 nm				
1-Q-1	~1	12.5	3.5	45	10	1	0.400	2	\$ 155.00
1-Q-2	2	12.5	4.5	45	10	2	0.700	2	155.00
1-Q-5	5	12.5	7.5	45	10	5	1.700	2	142.00
1-Q-10	10	12.5	12.5	45	10	10	3.500	2	92.00
1-Q-20	20	12.5	22.5	45	10	20	7.000	2	169.00
1-Q-40	40	12.5	42.5	45	10	40	14.00	2	228.00
1-Q-50	50	12.5	52.5	45	9.5	50	17.500	2	255.00
1-Q-100	100	12.5	102.5	45	9	100	35.000	2	412.00
Infrasil [®] or	equivalent Near In	ifrared Q	uartz wind	lows, Used	able range: 2	220 to 3800 n	m		
1-I-1	1	12.5	3.5	45	10	1	0.400	2	\$ 176.00
1-I-2	2	12.5	4.5	45	10	2	0.700	2	176.00
1-I-5	5	12.5	7.5	45	10	5	1.700	2	164.00
1-I-10	10	12.5	12.5	45	10	10	3.500	2	142.00
1-I-20	20	12.5	22.5	45	10	20	7.000	2	202.00
1-I-40	40	12.5	42.5	45	10	40	14.00	2	286.00
1-I-50	50	12.5	52.5	45	9.5	50	17.500	2	320.00
1-I-100	100	12.5	102.5	45	9	100	35.000	2	469.00

4


Standard Rectangular with stopper


Type 21 cells are identical to Standard Rectangular or Macro cells listed on page 4, except that instead of having an open rectangular hole at the top, a block is fused to the top of the cell with a ground hole to receive a PTFE stopper. This provides a more suitable seal for volatile liquids.

Cells with a path length of greater than 20mm are supplied with two stoppers and the height of the cell body is 48 mm with an overall height including the stopper of approximately 53mm.


Accessories:


Spacers for the 1, 2, and 5mm pathlength on page 7. Holders for Infrared instruments are available on page 39.

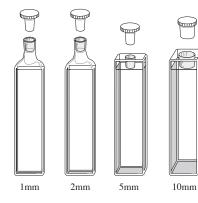

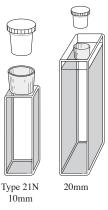

40mm

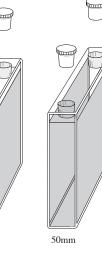

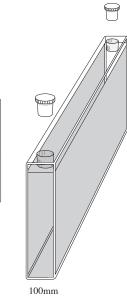

Catalog	Path Length		terior, n			or, mm	Nominal	Number	Polished	Price	
Number	mm	Width	Length	Height	Width	Length	Vol. ml	of Stoppers	Windows	per cell	
Optical Glass windows, Useable range: 334 to 2500 nm											
21-G-1	1	12.5	3.5	55	10	1	0.400	1	2	\$ 129.00	
21-G-2	2	12.5	4.5	55	10	2	0.700	1	2	129.00	
21-G-5	5	12.5	7.5	48	10	5	1.700	1	2	103.00	
21-G-10	10	12.5	12.5	48	10	10	3.500	1	2	90.00	
21-G-20	20	12.5	22.5	48	10	20	7.000	1	2	116.00	
21-G-40	40	12.5	42.5	48	10	40	14.00	2	2	150.00	
21-G-50	50	12.5	52.5	48	9.5	50	17.500	2	2	164.00	
21-G-100	100	12.5	102.5	48	9	100	35.000	2	2	312.00	
Special Optic	al Glass windows,	Useable r	ange: 32	0 to 2500 i	nm						
21-SOG-10	10	12.5	12.5	48	10	10	3.500	1	2	\$ 103.00	
Spectrosil [®] 1	Far UV Quartz win	idows, Use	eable ran	ge: 170 to	2700 nm						
21-Q-1	1	12.5	3.5	55	10	1	0.400	1	2	\$ 216.00	
21-Q-2	2	12.5	4.5	55	10	2	0.700	1	2	216.00	
21-Q-5	5	12.5	7.5	48	10	5	1.700	1	2	189.00	
21-Q-10	10	12.5	12.5	48	10	10	3.500	1	2	155.00	
21N-Q-10	10	12.5	12.5	44	10	10	3.000	1	2	189.00	
21-Q-20	20	12.5	22.5	48	10	20	7.000	1	2	242.00	
21-Q-40	40	12.5	42.5	48	10	40	14.00	2	2	325.00	
21-Q-50	50	12.5	52.5	48	9.5	50	17.500	2	2	372.00	
21-Q-100	100	12.5	102.5	48	9	100	35.000	2	2	652.00	
Infrasil [®] or	equivalent Near In	ifrared Qu	artz wind	dows, Use	able range: 2	220 to 3800	nm				
21-I-1	1	12.5	3.5	55	10	1	0.400	1	2	\$ 221.00	
21-I-2	2	12.5	4.5	55	10	2	0.700	1	2	221.00	
21-I-5	5	12.5	7.5	48	10	5	1.700	1	2	216.00	
21-I-10	10	12.5	12.5	48	10	10	3.500	1	2	202.00	
21-I-20	20	12.5	22.5	48	10	20	7.000	1	2	273.00	
21-I-40	40	12.5	42.5	48	10	40	14.00	2	2	404.00	
21-I-50	50	12.5	52.5	48	9.5	50	17.500	2	2	430.00	
21-I-100	100	12.5	102.5	48	9	100	35.000	2	2	692.00	


Semi-Micro

The purpose of Semi-Micro Type 9 cells is to reduce the sample volume necessary to make a measurement. The volume of a Semi-Micro cell is approximately 40% of an equivalent path length standard rectangular cell.

Semi-Micro cells are available with either clear walls parallel to the light beam or non transmitting self masking black walls which improve sensitivity by eliminating stray light from measurements. The interior sample compartment width of Semi-Micro cells is reduced to 4mm, therefore it is essential to determine that the light beam from the instrument passes centrally through the sample compartment. This is especially important for long path length cells. Both the Type 9 and Type 9B are supplied with a PTFE cover (not pictured with the Type 9B or 9B9).

Smaller sample volumes can be achieved by using cell Type 9/9 or 9B9 with a 9mm instead of a 3mm thick base. Note: Cells with a 9mm base are not suitable for an instrument with an 8.5mm 'Z' dimension.

Path Length

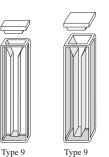

10mm

Type 9B

10mm


Interior, mm

5mm

Type 9B

5mm

Exterior, mm

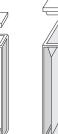


Type 9/9

10mm

Type 9B9

10mm

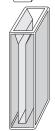


Type 9

20mm

Type 9B

20mm

Nominal

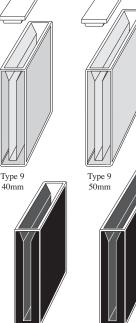

Type 9B

40mm

Polished



Type 9B

50mm

Price

Number Width Length Height Width Vol. ml Windows mm Length per cell Special Optical Glass windows, Useable range: 320 to 2500 nm **Clear Walls** 10 9-SOG-10 10 12.5 12.5 45 4 1.400 2 85.00 Spectrosil[®] Far UV Quartz windows, Useable range: 170 to 2700 nm 9-0-5 5 5 12.5 7.5 45 4 0.700 2 \$ 202.00 10 9-Q-10 12.5 45 10 2 12.5 4 1.400 164.00 9/9-Q-10 10 45 4 10 2 12.5 12.5 1.160 202.00 2 9-Q-20 20 12.5 22.5 45 4 20 2.800 247.00 40 4 2 9-Q-40 12.5 42.5 45 40 5.600 299.00 9-Q-50 50 12.5 52.5 45 4 50 7.000 2 325.00 Spectrosil[®] Far UV Quartz windows, Useable range: 170 to 2700 nm 9B-O-5 2 5 12.5 7.5 45 4 5 0.700 \$ 321.00 9B-Q-10 10 12.5 12.5 45 4 10 1.400 2 262.00 2 9B9-Q-10 10 12.5 12.5 45 4 10 1.160 335.00 9B-Q-20 20 12.5 22.5 45 4 20 2.800 2 397.00 9B-Q-40 40 12.5 42.5 45 4 40 5.600 2 523.00 9B-Q-50 50 12.5 52.5 45 4 50 7.000 2 652.00

Starna Cells, Inc (800) 228-4482



Catalog

Semi-Micro with stopper

Similar in all respects to the Type 9 Semi-Micro cells, with the exception that these cells are supplied with a PTFE stopper, or stoppers as pictured, to fit in a ground fitting to reduce evaporation of volatile samples (The Type 29B comes with stoppers even though the stoppers are not pictured). The Type 29 and Type 29B cells have an external height of 48 mm plus stopper.


Smaller sample volumes can be achieved by using cell Type 29/9 or 29B9 with a 9mm instead of a 3mm thick base.


Note: Cells with a 9mm base are not suitable for an instrument with an 8.5mm 'Z' dimension.

Short Pathlength Spacers

Black anodized aluminum spacers are available in three sizes for use with 1, 2, and 5mm path length cells, supporting them in a normal 12.5 x 12.5mm holder where there is no cell holder path length adjustment facility.

Cat. No Description Price, each for 1mm path length cells SPA-1 \$75.00 SPA-2 for 2mm path length cells 75.00 SPA-5 for 5mm path length cells 75.00

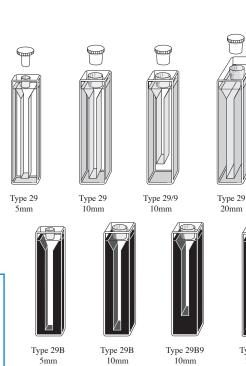


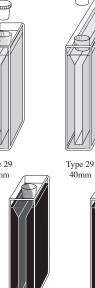

Type 29B

20mm

40mm

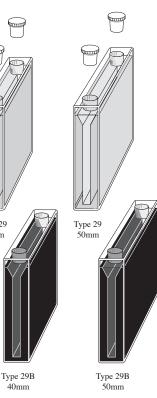

Catalog Number	Path Length		xterior, m Length		Interi Width	or, mm	Nominal Vol. ml	Polished Windows	Price
Number	mm	wiath	Length	neigni	wiam	Length	voi. mi	windows	per cell
Special Optic	al Glass windows,	Useable 1	range: 32	0 to 2500 i	nm			(lear Walls
29-SOG-10	10	12.5	12.5	48	4	10	1.400	2	150.00
Spectrosil [®] I	Far UV Quartz wir	idows, Us	eable ran	ge: 170 to	2700 nm				
29-Q-5	~ 5	12.5	7.5	48	4	5	0.700	2	\$242.00
29-Q-10	10	12.5	12.5	48	4	10	1.400	2	216.00
29/9-Q-10	10	12.5	12.5	48	4	10	1.160	2	273.00
29-Q-20	20	12.5	22.5	48	4	20	2.800	2	346.00
29-Q-40	40	12.5	42.5	48	4	40	5.600	2	412.00
29-Q-50	50	12.5	52.5	48	4	50	7.000	2	456.00
Spectrosil [®] I	Far UV Quartz wir	idows, Us	eable ran	ge: 170 to	2700 nm				
29B-Q-5	~ 5	12.5	7.5	48	4	5	0.700	2	\$ 423.00
29B-Q-10	10	12.5	12.5	48	4	10	1.400	2	343.00
29B9-Q-10	10	12.5	12.5	48	4	10	1.160	2	397.00
29B-Q-20	20	12.5	22.5	48	4	20	2.800	2	517.00
29B-Q-40	40	12.5	42.5	48	4	40	5.600	2	652.00
29B-Q-50	50	12.5	52.5	48	4	50	7.000	2	765.00


Micro Cells

Micro cells, Type 18, with a 2mm internal width are for use with sample volumes smaller than the standard rectangular or Semi-Micro cells. The nominal working volume of this range represents only 20% of the standard rectangular cells Type 1 with the same path length.

Micro cells are available with either clear walls or black walls which are opaque and improve sensitivity. The cells have an exterior height of 45mm. Each cell has an open top and is supplied with a PTFE cover (not pictured with the Type 18B or 18B9).

Smaller sample volumes can be achieved by using cell Type 18/9 or 18B9 with a 9mm instead of a 3mm thick base.

Note: Cells with a 9mm base are not suitable for instruments with an 8.5mm 'Z' dimension.

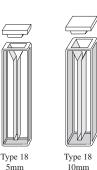

Type 18B

5mm




Type 18/9

10mm


Type 18B9

10mm

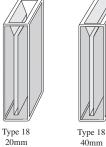

Type 18B

20mm

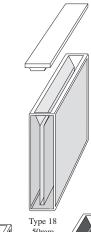

50mm

50mm



40mm



Catalog **Path Length** Exterior, mm Interior, mm Nominal Polished Price Number Width Length Height Width Length Vol. ml Windows per cell mm Special Optical Glass windows, Useable range: 320 to 2500 nm **Clear Walls** 18-SOG-10 10 12.5 45 2 10 0.700 2 98.00 12.5 Spectrosil[®] Far UV Quartz windows, Useable range: 170 to 2700 nm 18-Q-5 5 12.5 7.5 45 5 0.350 2 \$ 202.00 2 2 18-Q-10 10 12.5 12.5 45 2 10 0.700 176.00 2 18/9-Q-10 10 12.5 12.5 45 2 10 0.580 216.00 18-Q-20 20 12.5 22.5 45 2 20 1.400 2 233.00 40 12.5 42.5 45 2 40 2 18-Q-40 2.800 294.00 52.5 2 2 18-Q-50 50 12.5 45 50 3.500 320.00 Spectrosil[®] Far UV Quartz windows, Useable range: 170 to 2700 nm 18B-Q-5 5 12.5 7.5 45 5 0.350 2 \$356.00 2 18B-O-10 10 12.5 12.5 45 2 0.700 2 294.00 10 18B9-O-10 10 12.5 12.5 45 2 0.580 2 335.00 10 18B-Q-20 20 12.5 22.5 45 2 20 1.400 2 423.00 2 18B-Q-40 40 12.5 42.5 45 40 2.800 2 572.00 2 2 18B-Q-50 50 12.5 52.5 45 50 3.500 800.00

Type 18B

10mm

Micro Cells with stopper

Type 28 cells are identical to Type 18 Micro cells, except these cells are supplied with a PTFE stopper, or stoppers as pictured, to fit in a ground cone to reduce evaporation of volatile samples (The Type 28B comes with stoppers even though the stoppers are not pictured). The external height is 48mm plus the stopper. Micro cells are available with either clear walls which will transmit energy or opaque black walls for improved sensitivity. Alignment of Micro cells in instrument cell holders is critical.

Smaller sample volumes can be achieved by using cell Type 28/9 or 28B9 with a 9mm instead of a 3mm thick base.

Note: Cells with a 9mm base are not suitable for instruments with an 8.5mm 'Z' dimension.

Short Pathlength Spacers

Black anodized aluminum spacers are available in three sizes for use with 1, 2, and 5mm path length cells, supporting them in a normal 12.5 x 12.5mm holder where there is no cell holder path length adjustment facility. Cat. No Description Price, each

SPA-1for 1mm path length cells\$ 75.00SPA-2for 2mm path length cells75.00SPA-5for 5mm path length cells75.00


Type 28 5mm



Type 28B

10mm



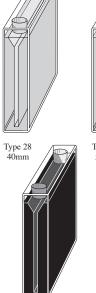

m



Type 28B 5mm

Type 28B9 10mm Type 28B 20mm Type 28B 40mm Type 28B 50mm

Catalog	Path Length	Ex	xterior, m	m	Inter	ior, mm	Nominal	Polished	Price			
Number	mm	Width	Length	Height	Width	Length	Vol. ml	Windows	per cell			
Spectrosil [®] H	Spectrosil [®] Far UV Quartz windows, Useable range: 170 to 2700 nm											
28-Q-5	5	12.5	7.5	48	2	5	0.350	2	\$ 281.00			
28-Q-10	10	12.5	12.5	48	2	10	0.700	2	228.00			
28/9-Q-10	10	12.5	12.5	48	2	10	0.580	2	260.00			
28-Q-20	20	12.5	22.5	48	2	20	1.400	2	346.00			
28-Q-40	40	12.5	42.5	48	2	40	2.800	2	430.00			
28-Q-50	50	12.5	52.5	48	2	50	3.500	2	464.00			
Spectrosil [®] H	Far UV Quartz win	dows, Use	eable ran	ge: 170 to	2700 nm							
28B-Q-5	5	12.5	7.5	48	2	5	0.350	2	\$ 477.00			
28B-Q-10	10	12.5	12.5	48	2	10	0.700	2	388.00			
28B9-Q-10	10	12.5	12.5	48	2	10	0.580	2	442.00			
28B-Q-20	20	12.5	22.5	48	2	20	1.400	2	590.00			
28B-Q-40	40	12.5	42.5	48	2	40	2.800	2	818.00			
28B-Q-50	50	12.5	52.5	48	2	50	3.500	2	894.00			

Sub-Micro Cells


Sub-Micro cells retain the exterior dimensions of a standard cell and are designed for the measurement of very small samples, with volumes ranging from $10\mu l$ to $160\mu l$. The entrance to the sample compartment is hemispherical, being designed without sharp corners to eliminate potential loss of sample due to capillary action.

The amount of sample required to fill the sample chamber is reduced to an absolute minimum, typically fifteen to twenty per cent greater than the absolute sample chamber volume.

The sample compartment in the Sub-Micro range of cells optimizes the alignment of the sample for maximum sensitivity in the spectrophotometer. The sample is easily inserted into and retrieved from the cell, by careful use of a pipette. The correct 'Z' dimension needs to be selected for this range, typically 8.5 or 15mm.

Cells are supplied with a PTFE cover and a tight sealing polyethylene cap.

'Z' Dimensio	on per Instrument
Manufacturer:	'Z' Dimension:
Agilent®	varies by instrument
Beckman [®]	8.5 mm
Bio-Rad [®]	8.5 mm
Eppendorf®	8.5 mm
GBC ^(R)	15 mm
Hewlett Packard	d [®] 15 mm
Hitachi®	varies by instrument
Jasco [®]	12 mm
Ocean Optics [®]	15 mm
Perkin-Elmer [®]	15 mm
Pharmacia [®]	15 mm
Shimadzu [®]	15 mm
StellarNet [®]	15 mm
Thermo Spectro	onic [®] 8.5 and 15 mm
Turner [®]	8.5 mm
Varian [®]	20 mm

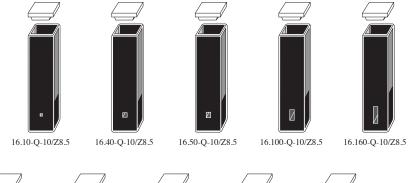

16.10-Q-10/Z15 16.40-Q-10/Z15 16.50-Q-10/Z15







16.100-Q-10/Z20





Catalog **Path Length** Exterior, mm Sample Chamber, mm 'Z' Dim. Nominal Polished Price Number Width Length Height Width Length Height Vol. ml Windows mm mm per cell Spectrosil[®] Far UV Quartz windows, Useable range: 170 to 2700 nm Sub-Micro, Open Top 16.10-Q-10/Z8.5 10 12.5 12.5 45 8.5 0.010 \$ 456.00 10 2 1 1 16.10-Q-10/Z15 10 2 10 12.5 12.5 45 15 0.010 456.00 1 1 16.10-Q-10/Z20 45 10 2 10 12.5 12.5 1 20 0.010 456.00 1 2 16.40-Q-10/Z8.5 10 2 10 2 0.040 12.5 12.5 45 8.5 364.00 16.40-Q-10/Z15 10 2 2 0.040 2 12.5 12.5 45 10 15 364.00 2 16.40-Q-10/Z20 10 12.5 12.5 45 2 10 2 20 0.040 364.00 2 16.50-Q-10/Z8.5 10 12.5 12.5 45 2 10 2.5 8.5 0.050 364.00 2 16.50-Q-10/Z15 10 12.5 12.5 45 2 10 2.5 15 0.050 364.00 16.50-O-10/Z20 2 2 10 12.5 12.5 45 10 2.5 20 0.050 364.00 2 2 16.100-Q-10/Z8.5 10 12.5 12.5 45 10 5 8.5 0.100 364.00 2 2 16.100-Q-10/Z15 10 12.5 12.5 45 10 5 15 0.100 364.00 12.5 2 2 16.100-Q-10/Z20 12.5 45 5 364.00 10 10 20 0.100 2 2 16.160-Q-10/Z8.5 12.5 12.5 45 10 8 8.5 364.00 10 0.160 2 8 2 16.160-Q-10/Z15 10 12.5 12.5 45 10 15 0.160 364.00 2 16.160-Q-10/Z20 10 12.5 12.5 45 2 10 8 20 0.160 364.00

10

Low Head Space Cells, Thermal Melt

Low head space Sub-Micro cells are designed for thermal melt. The low head space cells hold volumes ranging from 10μ l to 160μ l while retaining standard exterior dimensions. The entrance to the sample compartment is hemispherical, being designed without sharp corners to eliminate potential loss of sample due to capillary action.


Cells are supplied with a special PTFE stopper. The liquid tight stopper reduces the volume of air by 95% over the volume of the sub-micro cells on the previous page.

The sample compartment in the Sub-Micro range of cells optimizes the alignment of the sample for maximum sensitivity in the spectrophotometer. The sample is easily inserted into and retrieved from the cell, by careful use of a pipette. The correct 'Z' dimension needs to be selected for this range, typically 8.5, 15 or 20mm.

The amount of sample required to fill the sample chamber is reduced to an absolute minimum, typically fifteen to twenty per cent greater than the absolute sample chamber volume.



26.160LHS-Q-10/Z15

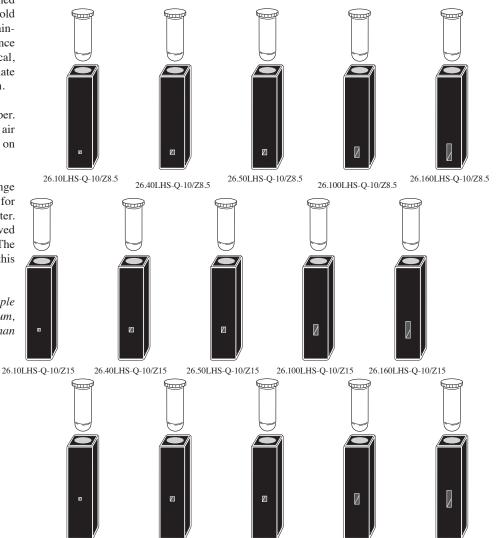

26.10LHS-Q-10/Z20

26.40LHS-Q-10/Z20 26.50LHS-Q-10/Z20

26.100LHS-Q-10/Z20


Catalog I Number	Path Lengt mm		xterior, r Length		-	e Chambe Length	·	'Z' Dim mm	Nominal Vol. ml	Polished Windows	Price per cell
Spectrosil [®] Far U	UV Quartz	windows	, Useabl	e range: i	170 to 27	00 nm			Ult	ra-Micro, Pi _l	oette filling
26.10LHS-Q-10/Z8	3.5 10	12.5	12.5	45	1	10	1	8.5	0.010	2	\$ 706.00
26.10LHS-Q-10/Z1	5 10	12.5	12.5	45	1	10	1	15	0.010	2	706.00
26.10LHS-Q-10/Z2	0 10	12.5	12.5	45	1	10	1	20	0.010	2	706.00
26.40LHS-Q-10/Z8	5 10	12.5	12.5	45	2	10	2	8.5	0.040	2	600.00
26.40LHS-Q-10/Z1	5 10	12.5	12.5	45	2	10	2	15	0.040	2	600.00
26.40LHS-Q-10/Z2	2 0 10	12.5	12.5	45	2	10	2	20	0.040	2	600.00
26.50LHS-Q-10/Z8	5 10	12.5	12.5	45	2	10	2.5	8.5	0.050	2	600.00
26.50LHS-Q-10/Z1	5 10	12.5	12.5	45	2	10	2.5	15	0.050	2	600.00
26.50LHS-Q-10/Z2	2 0 10	12.5	12.5	45	2	10	2.5	20	0.050	2	600.00
26.100LHS-Q-10/Z	8.5 10	12.5	12.5	45	2	10	5	8.5	0.100	2	600.00
26.100LHS-Q-10/Z	2 15 10	12.5	12.5	45	2	10	5	15	0.100	2	600.00
26.100LHS-Q-10/Z	2 0 10	12.5	12.5	45	2	10	5	20	0.100	2	600.00
26.160LHS-Q-10/Z	8.5 10	12.5	12.5	45	2	10	8	8.5	0.160	2	600.00
26.160LHS-Q-10/Z	15 10	12.5	12.5	45	2	10	8	15	0.160	2	600.00
26.160LHS-Q-10/Z	2 0 10	12.5	12.5	45	2	10	8	20	0.160	2	600.00

Screw Cap Cells, Closed and Septum - GL14

Cell types with suffix **GL14** are able to be used under anaerobic conditions. Cells designated with **GL14** have a standard GL14 threaded top and cap which can be either a septum cap GL14-**S** or a plain closed cap GL14-**C**. Each cap contains a silicon seal to give you an air tight fit. The external diameter of the cap along the bottom rim is close to 20mm.




Type 1 10mm GL14-C

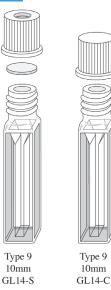



10mm

GL14-S

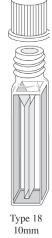

GL14-C

GL14 Caps, Closed & Septum

Additional screw caps to fit GL14 threaded cells are available as either the closed cap or septum seal cap.

Cat. No	Description	Price, each
GL14-C	Closed cap	\$ 6.00
GL14-S	Septum seal cap	6.00
GL14/SI	Septum Insert	3.00



Catalog	P	ath Length	n E	xterior, 1	nm	Interi	or, mm	Nominal	Polished	Price
Number	Description	mm	Width	Length	Height	Width	Length	Vol. ml	Windows	per cell
Screw Cap Cells	with Closed Cap:									
Special Optical Glas	ss windows, Useable ra	ange: 320 to	o 2500	nm						
1-SOG-10-GL14-C	Closed Cap, Standard	10	12.5	12.5	68	10	10	3.500	2	\$ 202.00
Spectrosil [®] Far UV	Quartz windows, Used	able range:	170 to	2700 nn	ı					
1-Q-10-GL14-C	Closed Cap, Standard	10	12.5	12.5	68	10	10	3.500	2	\$ 228.00
9-Q-10-GL14-C	Closed Cap, Semi-Mic	ro 10	12.5	12.5	68	4	10	1.400	2	320.00
18-Q-10-GL14-C	Micro	10	12.5	12.5	68	2	10	0.700	2	359.00
Screw Cap Cells	with Septum Caps	•								
Special Optical Glas	ss windows, Useable ra	ange: 320 to	o 2500	nm						
1-SOG-10-GL14-S	Septum Cap, Standard	10	12.5	12.5	68	10	10	3.500	2	\$ 202.00
Spectrosil [®] Far UV	Quartz windows, Used	able range:	170 to	2700 nn	ı					
1-Q-10-GL14-S	Septum Cap, Standard	10	12.5	12.5	68	10	10	3.500	2	\$ 228.00
9-Q-10-GL14-S	Septum Cap, Semi-Mic	cro 10	12.5	12.5	68	4	10	1.400	2	320.00
18-Q-10-GL14-S	Septum Cap, Micro	10	12.5	12.5	68	2	10	0.700	2	359.00

Screw Cap Cells, Closed and Septum - GL14

Cell types with suffix **GL14** are able to be used under anaerobic conditions. Cells designated with **GL14** have a standard GL14 threaded top and cap which can be either a septum cap GL14-**S** or a plain closed cap GL14-**C**. Each cap contains a silicon seal to give you an air tight fit.

The black walled versions of these cells provide improved sensitivity as the cell insures that all the light that strikes your instrument's photo detector has traveled through your sample.


	Screw caps to either closed of Cat. No GL14-C GL14-S	cap or septum sea Description Closed cap Septum sea	ed cells. Available al cap. n Price, eac \$6.0 l cap 6.0	h 0 0					
	GL14/SI	Septum Inse	ert 3.0	0	Type 9B 10mm GL14-S	Type 9 10mi GL14	m 10	0mm	Гуре 18 10mm GL14-C
Catalog Number	Description	Path Length mm	n Exterior, n Width Length			or, mm Length	Nominal Vol. ml	Polished Windows	
	ls with Closed Ca	p:	8	8		0			L
S <i>pectrosil[®] Far U</i> 28-Q-10-GL14-C 28B-Q-10-GL14-C	V Quartz windows, U Semi-Micro, Black Micro, Black Wall		: 170 to 2700 nm 12.5 12.5 12.5 12.5	68 68	4 2	10 10	1.400 0.700	2 2	437.00 491.00
	Is with Septum Co V Quartz windows, U Semi-Micro, Black Micro, Black Wall	Jseable range:	: 170 to 2700 nm 12.5 12.5 12.5 12.5	68 68	4 2	10 10	1.400 0.700	2 2	437.00 491.00
	p Cells, Close	-		-					
Cell types with suffi small ST cap have t cap, ST-C. Each c	p Cells, Close ix ST are able to be use hreaded top and cap whi ap contains a silicon sea end beyond the perimete	ed under anaerol ich can be either al to give you	bic conditions. Cel a septum cap, ST-	ls desigi S or a pl	nated with ain closed				
Cell types with suffi small ST cap have to cap, ST-C. Each co diameter) do not extra Small Scr Additonal scr as either the co Cat. No ST-C Co ST-S S	ix ST are able to be use hreaded top and cap whi ap contains a silicon sea	ed under anaerol ich can be either al to give you a er of the cell. & Septum ed cells are avail cap. Price, each \$ 29.00 34.00	bic conditions. Cel a septum cap, ST- an air tight fit. Th	ls design S or a pl le ST ca	nated with ain closed		Type 10mr ST-S	m	Type 1 10mm ST-C
Cell types with suffi small ST cap have to cap, ST-C. Each co diameter) do not ext Small Scr Additonal scr as either the co Cat. No ST-C C ST-S S ST-SI S Catalog	ix ST are able to be use hreaded top and cap whi ap contains a silicon sea end beyond the perimete rew Caps, Closed ew caps to fit ST threade closed cap or septum seal Description Closed cap Septum seal cap	ed under anaerol ich can be either al to give you a er of the cell. & Septum ed cells are avail cap. Price, each \$ 29.00 34.00	bic conditions. Cel a septum cap, ST- an air tight fit. Th lable	ls design S or a pl te ST ca	nated with ain closed ps (12mm ST-C mall Solid Cap	Dr, mm	10mr ST-S	m	10mm ST-C Price
Cell types with suffi small ST cap have to cap, ST-C. Each co diameter) do not extr Small Scr Additonal scr as either the co Cat. No ST-C Co ST-S ST-SI ST-SI ST-SI Catalog Number Screw Cap Cell	ix ST are able to be use hreaded top and cap whi ap contains a silicon sea end beyond the perimete rew Caps, Closed ew caps to fit ST threade closed cap or septum seal Description Closed cap Septum seal cap Septum Insert, pack of 10	ed under anaerol ich can be either al to give you a er of the cell. & Septum ed cells are avail cap. Price, each \$ 29.00 34.00 0 25.00 Path Length mm	bic conditions. Cel a septum cap, ST- an air tight fit. Th lable ST-S Small Septum Ca Exterior, n Width Length	ls design S or a pl le ST ca p S m Height	nated with ain closed ps (12mm ST-C mall Solid Cap	Dr, mm	10mr ST-S	m S Polished	10mm ST-C Price
Cell types with suffi small ST cap have to cap, ST-C. Each co diameter) do not extr Additonal scr as either the co Cat. No ST-C C ST-S S ST-SI S Catalog Number Screw Cap Cell Spectrosil [®] Far U I-Q-10-ST-C	ix ST are able to be use hreaded top and cap whi ap contains a silicon sea end beyond the perimeter rew Caps, Closed rew caps to fit ST threade closed cap or septum seal Description Closed cap Septum seal cap Septum Insert, pack of 10 Description Is with Closed Cap W Quartz windows, U	ed under anaerol ich can be either al to give you a er of the cell. & Septum ed cells are avail cap. Price, each \$ 29.00 34.00 0 25.00 Path Length mm p: Useable range: Useable range:	bic conditions. Cel a septum cap, ST- an air tight fit. Th lable ST-S Small Septum Ca Exterior, n Width Length 1 : 170 to 2700 nm 12.5 12.5	ls design S or a pl le ST ca be ST ca s m Height 58	nated with ain closed ps (12mm ST-C mall Solid Cap Interio Width	or, mm Length	IOmr ST-S Nominal Vol. ml	Polished Windows	ST-C Price per cel

13

Short Path Length, Demountable

The Type 20 cell is designed for use in applications where a short path length and large window surface area is needed. This cell design is especially useful for opaque or viscose solutions. The cell consists of one fully fused component combining one window and the sides of the cell which forms the path length. The other window is separate and the cell can therefore be taken apart for cleaning. The cell wall and the removable window are polished so flat that the two will form a seal when placed in contact. Viscous and aqueous samples will allow the best seal to be formed, organic solvents will creep into the seal area for a less reliable seal.

Type 20/O is open at one end similar to a Standard cell, unlike the Type 20/C which is completely enclosed when the window is in place. It is recommended that the cell holder Type CH-2049, shown below, is used with these cells to keep the cell together while standing in the instrument cell holder.

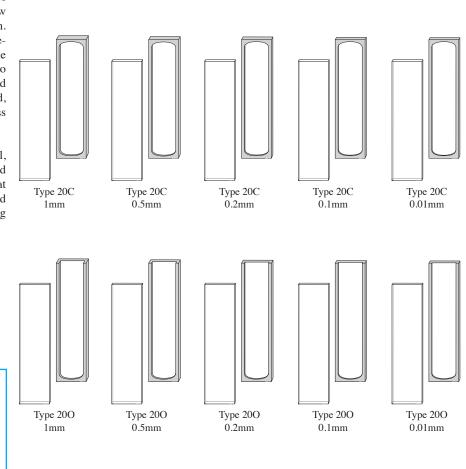

Cell Holder

This cell holder is designed to lightly clamp the type 20 cells in a vertical orientation for use in a standard instrument cell holder. Exterior dimensions are 12.5x 12.5 x 52mm. The base is 9mm thick. This cell holder works with cells having a 2mm pathlength or less.



Cat. NoDescriptionPrice, eachCH-2049Cell holder\$ 103.00

Catalog Number	Path Length mm		terior, n Length		1	Chamb Length	/	Nominal Vol ml	Polished Windows	Price per cell
Spectrosil [®]	Far UV Quartz v	vindows	, Useable	range: 1	70 to 2700 nn	ı				Open Top
20/O-Q-0.01	0.01	12.5	2.5	45	10	0.01	40	0.004	2	\$ 372.00
20/O-Q-0.1	0.1	12.5	2.6	45	10	0.1	40	0.040	2	247.00
20/O-Q-0.2	0.2	12.5	2.7	45	10	0.2	40	0.080	2	247.00
20/O-Q-0.5	0.5	12.5	3	45	10	0.5	40	0.190	2	247.00
20/O-Q-1	1	12.5	3.5	45	10	1	40	0.390	2	228.00
									Closed	Both Ends
20/C-Q-0.01	0.01	12.5	2.5	45	8	0.01	40	0.003	2	\$ 372.00
20/C-Q-0.1	0.1	12.5	2.6	45	8	0.1	40	0.030	2	255.00
20/C-Q-0.2	0.2	12.5	2.7	45	8	0.2	40	0.060	2	255.00
20/C-Q-0.5	0.5	12.5	3	45	8	0.5	40	0.150	2	255.00
20/C-Q-1	1	12.5	3.5	45	8	1	40	0.310	2	228.00

14

Short Path Length, Flow Cells

Cell Types 48 and 49 are designed for use in various applications where a flow cell with a small path length and a large window surface area is needed. Flow tubes are fused to the front window of the cell, in the case of Type 49, the removable window. The cells will fit any spectrophotometer which can accommodate a normal cell with a 40mm path length. This allows introduction of the inlet/outlet tubes parallel to the light beam. The Type 48 is a fully fused cell.

The demountable cells Type 49, with removable window, must be held together when being used. Cell holder Type CH-2049, shown below, is designed for this purpose. These cells may be taken apart for easy cleaning. This cell design is especially useful for opaque or viscose solutions. The cell consists of one fully fused component combining one window and the sides of the cell which forms the path length. The other window is separate and the cell can therefore be taken apart for cleaning. The cell wall and the removable window are polished so flat that the two will form a seal when placed in contact. Viscous and aqueous samples will allow the best seal to be formed, organic solvents will creep into the seal area for a less reliable seal.

Cell Holder

This cell holder is designed to lightly clamp the type 20 cells in a vertical orientation for use in a standard instrument cell holder. Exterior dimensions are 12.5 x 12.5 x 52mm. The base is 9mm thick. This cell holder works with cells having a 2mm pathlength or less.

Cat. No **CH-2049**



Description Price, each \$ 103.00 Cell holder

Type 48 Side View X-Section Type 48 Type 48 Type 48 Type 48 Type 48 2mm 5mm 1mm 0.5mm 0.2mm 6 6 Type 49 Side View X-Section Type 49 Type 49 Type 49 Type 49 5mm 2mm 1mm

Flow cell Tubul	ation dimensi	ions:
OD	ID	Length
4.2mm	2.0mm	16mm

Catalog Number	Path Length mm		terior, 1 Length	nm Height	-	e Chaml Length	ber, mm Height	Volume in ml	Polished Windows	Price per cell
Spectrosil [®] H	Far UV Quartz v	vindows	, Useabl	e range: 1	70 to 2700 ni	n				
48-Q-0.1	0.1	12.5	2.6	45	8	0.1	40	0.030	2	\$ 351.00
48-Q-0.2	0.2	12.5	2.7	45	8	0.2	40	0.060	2	351.00
48-Q-0.5	0.5	12.5	3	45	8	0.5	40	0.150	2	351.00
48-Q-1	1	12.5	3.5	45	8	1	40	0.300	2	339.00
48-Q-2	2	12.5	4.5	45	8	2	40	0.600	2	339.00
48-Q-5	5	12.5	7.5	45	8	5	40	1.560	2	351.00
Spectrosil [®] F	Far UV Quartz v	vindows	, Useabl	e range: 1	70 to 2700 ni	n				
49-Q-0.1	0.1	12.5	2.6	45	8	0.1	40	0.030	2	\$ 351.00
49-Q-0.2	0.2	12.5	2.7	45	8	0.2	40	0.060	2	351.00
49-Q-0.5	0.5	12.5	3	45	8	0.5	40	0.150	2	351.00
49-Q-1	1	12.5	3.5	45	8	1	40	0.300	2	339.00
49-Q-2	2	12.5	4.5	45	8	2	40	0.600	2	339.00
49-Q-5	5	12.5	7.5	45	8	5	40	1.560	2	351.00

www.starnacells.com


Type 48

0.1mm

0.5mm

Flow Cells, In-Line or Microscope Analysis

The type 45 flow cells provide a large surface to view materials as they flow through the chamber. The thin dimension (3 to 4.5mm external thickness) makes this cell ideal for use with most microscopes.

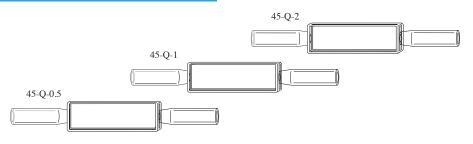

The inlet/outlet tube dimensions are 25mm long: Up to 2mm pathlength: 2mm ID, 4mm OD 5mm pathlength: 5mm ID, 7mm OD


Catalog Number	Path Length mm		terior, Lengtl	mm h Height	Interi Width	ior, mm Length	Nominal Vol. ml	Polished Windows	Price per cell
Spectrosil [®] 1	Far UV Quartz win	dows, Use	able ra	inge: 170 to 2	2700 nm				
45-Q-0.5	0.5	12.5	3.0	40 + tubes	10	0.5	0.225	2	\$482.00
45-Q-1	1	12.5	3.5	40 + tubes	10	1	0.450	2	398.00
45-Q-2	2	12.5	4.5	40 + tubes	10	2	0.900	2	372.00

Flow Cells, Standard & Semi-Micro

Type 71 Standard, and 72 Semi-Microseries flow cells are suitable for use with instruments which have either 8.5 or 15 mm 'Z' dimensions, where the volume of the sample available is not a critical factor.

Type 46 Standard and 47 Semi-Micro are for in-line continuous flow type applications. Samples normally flow from the bottom to the top.

Flow cell Tubulation dimensions:									
OD	ID	Length							
3.5mm	2.2mm	10mm							

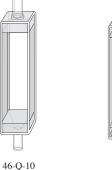

10

1

10

12.5

12.5

12.5

12.5

12.5

12.5

48


48


48


7


4

4

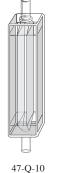

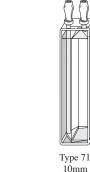

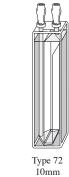

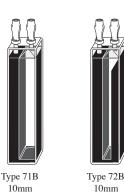

10mm

3.000

0.180

1.800

2


2

2

577.00

786.00

585.00

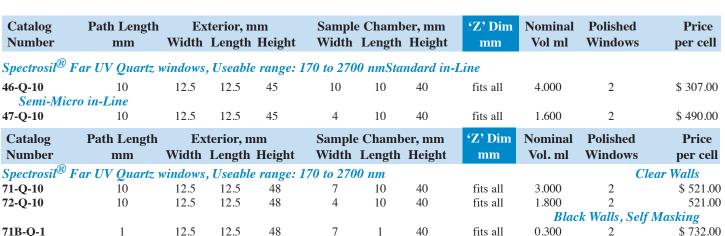

40

40

40

10

1

10

10mm

Catalog

Number

46-O-10

47-Q-10

Catalog

Number

71-O-10

72-Q-10

71B-O-1

71B-Q-10

72B-Q-1

72B-Q-10

Starna Cells, Inc (800) 228-4482

fits all


fits all

fits all

M6 Threaded Flow Cells

The fully fused Far UV quartz flow cell body mounted in a threaded anodized aluminum jacket. The aluminum exterior has internally threaded (M6 threading) inlet/outlet ports into which the PTFE flow tubes are connected. This allows for a tight seal between the inlet/outlet tubes and the flow cell body while making connections very easy to change. All cells and tubing connections are tested to 5 bar pressure (75 psi).

The exterior dimensions of the flow cells are made to international standards which should allow their use with commonly available commercial instruments. Almost any UV/Vis spectrophotometer has a cell holder which will work with the flow cells having a pathlength of 10mm or less. For pathlengths greater than 10mm consult your instrument handbook for the correct cell holder to use from your instrument manufacturer.

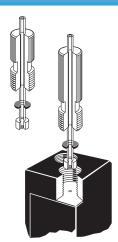

Cut away showing M6 threaded connectors


M6 Flow cell Tubulation dimensions: *OD: 1.6mm ID:1mm Length: 400mm*

Flow Cells, Medium Aperture

The Type 583.4 series of flow cells are designed to fit most spectrophotometers and are available in most 'Z' dimensions. They have a rectangular window in a black self masking configuration and are widely used for many standard applications including Tablet dissolution. The sample chamber width is 4mm and is available with path lengths of 5, 10, 20, 40, 50, and 100mm. The inlet and outlet tubes are connected into the flow cell with a M6 threaded connector. The PTFE tubing has an outer diameter of about 1.6mm and an inner diameter of about 1mm. The length of the tubing is approximately 400mm. Each cell is supplied with an inlet and an outlet tube. Extra tubing sets may be purchased.

'Z' Dimension per Instrument

Manufacturer:	'Z' Dimension:
Agilent®	varies by instrument
Beckman [®]	8.5 mm
Bio-Rad [®]	8.5 mm
Eppendorf [®]	8.5 mm
GBC ^(R)	15 mm
Hewlett Packard	15 mm
Hitachi [®]	varies by instrument
Jasco [®]	12 mm
Ocean Optics [®]	15 mm
Perkin-Elmer [®]	15 mm
Pharmacia®	15 mm
Shimadzu®	15 mm
StellarNet [®]	15 mm
Thermo Spectron	ic [®] 8.5 and 15 mm
Turner [®]	8.5 mm
Varian [®]	20 mm

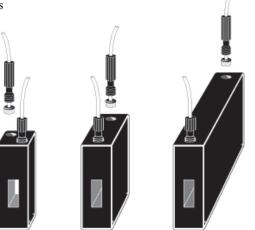


Catalog Number	Path Length mm		terior, n Length	ım Height	-	le Chaml Length		'Z' Dim mm	Nominal Vol. ml	Polished Windows	Price per cell
Spectrosil [®] Far	UV Quartz wi	indow, U	seable ra	inge: 170 i	to 2700 r	ım					
583.4-Q-5/Z8.5	5	12.5	12.5	35	4	5	11	8.5	0.225	2	\$ 678.00
583.4-Q-10/Z8.	5 10	12.5	12.5	35	4	10	11	8.5	0.450	2	516.00
583.4-Q-20/Z8.4	5 20	12.5	22.5	35	4	20	11	8.5	0.900	2	874.00
583.4-Q-40/Z8.4	5 40	12.5	42.5	35	4	40	11	8.5	1.800	2	1,078.00
583.4-Q-50/Z8.	5 50	12.5	52.5	35	4	50	11	8.5	2.250	2	1,240.00
583.4-Q-100/Z8	3.5 100	12.5	102.5	35	4	100	11	8.5	4.500	2	1,789.00
583.4-Q-5/Z15	5	12.5	12.5	35	4	5	11	15	0.225	2	\$ 678.00
583.4-Q-10/Z15	10	12.5	12.5	35	4	10	11	15	0.450	2	516.00
583.4-Q-20/Z15	20	12.5	22.5	35	4	20	11	15	0.900	2	874.00
583.4-Q-40/Z15	40	12.5	42.5	35	4	40	11	15	1.800	2	1,078.00
583.4-Q-50/Z15	50	12.5	52.5	35	4	50	11	15	2.250	2	1,240.00
583.4-Q-100/Z1	5 100	12.5	102.5	35	4	100	11	15	4.500	2	1,789.00
583.4-Q-5/Z20	5	12.5	12.5	35	4	5	11	20	0.225	2	\$ 678.00
583.4-Q-10/Z20	10	12.5	12.5	35	4	10	11	20	0.450	2	516.00
583.4-Q-20/Z20	20	12.5	22.5	35	4	20	11	20	0.900	2	874.00
583.4-Q-40/Z20	40	12.5	42.5	35	4	40	11	20	1.800	2	1,078.00
583.4-Q-50/Z20	50	12.5	52.5	35	4	50	11	20	2.250	2	1,240.00
583.4-Q-100/Z2	2 0 100	12.5	102.5	35	4	100	11	20	4.500	2	1,789.00
M6-SET	Set of e	xtra inlet/	outlet tub	es with M6	connecto	ors (two tu	bes per set)				111.00





Flow Cells, Wide Aperture

The Type 583.65 series of flow cells have a wide aperture (6.5mm) that are optimized for instruments with a wide light beam. They have a rectangular window in a black self masking configuration and are widely used for many standard applications including Tablet dissolution. The sample chamber width of 6.5mm is available with path lengths of 0.1, 0.2, 0.5, 1, 2, 5, 10, 20, 40, 50, and 100mm. The inlet and outlet tubes are connected into the flow cell with a M6 threaded connector. The PTFE tubing has an outer diameter of about 1.6mm and an inner diameter of about 1mm. The length of the tubing is approximately 400mm. Each cell is supplied with an inlet and an outlet tube. Extra tubing sets may be purchased.

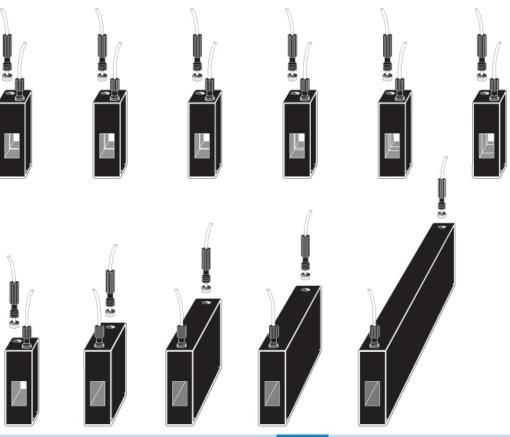

Catalog Number	Path Length mm		terior, n Length		-	e Chamb Length		'Z' Dim mm	Nominal Vol. ml	Polished Windows	Price per cell
Spectrosil [®] Fa	r UV Ouartz,	window,	Useable	range: 12	70 to 270) nm					
583.65-Q-0.1/Z1		12.5	12.5	35	6.5	0.1	11	15	0.029	2	\$ 986.00
583.65-Q-0.2/Z1		12.5	12.5	35	6.5	0.2	11	15	0.036	2	986.00
583.65-Q-0.5/Z1	5 0.5	12.5	12.5	35	6.5	0.5	11	15	0.072	2	92600
583.65Q-1/Z15	1	12.5	12.5	35	6.5	1	11	15	0.072	2	725.00
583.65-Q-2/Z15	2	12.5	12.5	35	6.5	2	11	15	0.290	2	712.00
583.65-Q-5/Z15	5	12.5	12.5	35	6.5	5	11	15	0.360	2	634.00
583.65-Q-10/Z15	10	12.5	12.5	35	6.5	10	11	15	0.720	2	542.00
583.65-Q-20/Z15	20	12.5	22.5	35	6.5	20	11	15	1.400	2	817.00
583.65-Q-40/Z15		12.5	42.5	35	6.5	40	11	15	2.900	2	1,157.00
583.65-Q-50/Z15	50	12.5	52.5	35	6.5	50	11	15	3.600	2	1,254.00
583.65-Q-100/Z1	5 100	12.5	102.5	35	6.5	100	11	15	7.200	2	1,679.00
583.65-Q-0.1/Z2	0 0.1	12.5	12.5	40	6.5	0.1	11	20	0.029	2	\$ 986.00
583.65-Q-0.2/Z2		12.5	12.5	40	6.5	0.2	11	20	0.036	2	986.00
583.65-Q-0.5/Z2	0 0.5	12.5	12.5	40	6.5	0.5	11	20	0.072	2	926.00
583.65Q-1/Z20		12.5	12.5	40	6.5	1	11	20	0.072	2	725.00
583.65-Q-2/Z20	2	12.5	12.5	40	6.5	2	11	20	0.290	2	712.00
583.65-Q-5/Z20	5	12.5	12.5	40	6.5	5	11	20	0.360	2	634.00
583.65-Q-10/Z20		12.5	12.5	40	6.5	10	11	20	0.720	2	542.00
583.65-Q-20/Z20	20	12.5	22.5	40	6.5	20	11	20	1.400	2	817.00
583.65-Q-40/Z20		12.5	42.5	40	6.5	40	11	20	2.900	2	1,157.00
583.65-Q-50/Z20		12.5	52.5	40	6.5	50	11	20	3.600	2	1,254.00
583.65-Q-100/Z2		12.5	102.5	40	6.5	100	11	20	7.200	2	1,679.00
M6-SET	Set of e	xtra inlet/	outlet tub	es with Mo	6 connecto	rs (two tul	bes per set)				111.00



Flow Cells, Wide Aperture, short window

The Type 583.65.65 series of flow cells have a wide aperture (6.5mm) and a smaller sample size. They have a rectangular window in a black self masking configuration and are widely used for many standard applications including Tablet dissolution. The sample chamber width of 6.5mm is available with path lengths of 1, 2, 5, 10mm. The inlet and outlet tubes are connected into the flow cell with a M6 threaded connector. The PTFE tubing has an outer diameter of about 1.6mm and an inner diameter of about 1mm. The length of the tubing is approximately 400mm. Each cell is supplied with an inlet and an outlet tube. Extra tubing sets may be purchased.



Catalog I Number	Path Length mm		terior, n Length	nm Height	1		ber, mm Height	'Z' Dim mm	Nominal Vol. ml	Polished Windows	Price per cell
Spectrosil [®] Fa	r UV Quartz	window,	Useable	range: 1	70 to 2700) nm					
583.65.65-Q-1/Z	15 1	12.5	12.5	35	6.5	1	6.5	15	0.042	2	\$ 744.00
583.65.65-Q-2/Z	15 2	12.5	12.5	35	6.5	2	6.5	15	0.084	2	731.00
583.65.65-Q-5/Z	15 5	12.5	12.5	35	6.5	5	6.5	15	0.210	2	652.00
583.65.65-Q-10/2	Z15 10	12.5	12.5	35	6.5	10	6.5	15	0.420	2	568.00
583.65.65-Q-1/Z	20 1	12.5	12.5	40	6.5	1	6.5	20	0.042	2	\$ 769.00
583.65.65-Q-2/Z	20 2	12.5	12.5	40	6.5	2	6.5	20	0.084	2	756.00
583.65.65-Q-5/Z	20 5	12.5	12.5	40	6.5	5	6.5	20	0.210	2	678.00
583.65.65-Q-10/2	Z20 10	12.5	12.5	40	6.5	10	6.5	20	0.420	2	587.00
M6-SET	Set of e	xtra inlet	outlet tub	bes with M	6 connector	rs (two tu	ibes per set)				111.00

Flow Cells, Short Path Length

This series of flow cell is typified by a tall sample compartment with a window height of 17.5mm for use in both 8.5 and 15mm 'Z' dimension situations. Path lengths as short as 0.01mm are available in this range. Smooth flow characteristics at very short path lengths are achieved with a bypass facility at each side of the window which also avoids increased pressure in the cell and helps to minimize air bubbles. The width of the sample chamber is 4mm. The PTFE inlet and outlet tubes are connected into the flow cell with a M6 threaded connector. The tubing has an outer diameter of about 1.6mm and an inner diameter of about 1mm. The length of the tubing is approximately 400mm. Each cell is supplied with an inlet and an outlet tube. Extra tubing sets may be purchased.




19

Catalog Number	Path Length mm		xterior, r Length	nm Height	-		ber, mm 1 Height	'Z' Dim mm	Nominal Vol. ml	Polished Windows	Price per cell
Spectrosil [®] Fe	ar UV Quartz w	indow, l	Useable i	range: 170	to 2700	nm					
584.4-Q-0.01	0.01	12.5	12.5	35	4	0.01	17.5	8.5 & 15	0.036	2	\$ 1,104.00
584.4-Q-0.05	0.05	12.5	12.5	35	4	0.05	17.5	8.5 & 15	0.039	2	1,044.00
584.4-Q-0.1	0.1	12.5	12.5	35	4	0.1	17.5	8.5 & 15	0.041	2	882.00
584.4-Q-0.2	0.2	12.5	12.5	35	4	0.2	17.5	8.5 & 15	0.047	2	882.00
584.4-Q-0.5	0.5	12.5	12.5	35	4	0.5	17.5	8.5 & 15	0.095	2	882.00
584.4-Q-1	1	12.5	12.5	35	4	1	17.5	8.5 & 15	0.120	2	686.00
584.4-Q-2	2	12.5	12.5	35	4	2	17.5	8.5 & 15	0.240	2	686.00
M6-SET	Set of ex	tra inlet/	outlet tube	es with M6	connectors	s (two tu	bes per set)				111.00

Flow Cells, Ultra-micro, round aperture

Type 585.1 and 585.15 are designed to optimize a very small sample volume by using a cylindrical sample chamber with a very small diameter. The diameter of the aperture on the 585.1 series of flow cells is 1 mm and the aperture on the 585.15 series is 1.5mm. The choice of the correct 'Z' dimension for your instrument is critical as you must ensure that the light beam of the instrument passes through the sample chamber of the cell. The cell is designed with a threaded jacket which will accept M6 threaded connectors. Available 'Z' Dimensions are 8.5, 15 (cell height 35mm) and 20mm (cell height 40mm).

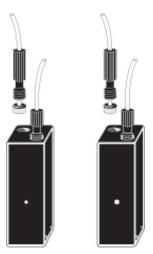

Catalog Number	Path Length mm		terior, n Length	nm Height	Sample Chan Diameter Ø	/	'Z' Dim mm	Nominal Vol. ml	Polished Windows	Price per cell
Spectrosil [®] Fa	r UV Quartz, w	vindows,	Useable	e range: 17	70 to 2700 nm				M6 Thread	ed Jacket
585.1-Q-10/Z8.5	5 10	12.5	12.5	35	1	10	8.5	0.008	2	\$ 626.00
585.15-Q-10/Z8	.5 10	12.5	12.5	35	1.5	10	8.5	0.018	2	573.00
585.1-Q-10/Z15	10	12.5	12.5	35	1	10	15	0.008	2	626.00
585.15-Q-10/Z1	5 10	12.5	12.5	35	1.5	10	15	0.018	2	573.00
585.1-Q-10/Z20	10	12.5	12.5	40	1	10	20	0.008	2	665.00
585.15-Q-10/Z2	0 10	12.5	12.5	40	1.5	10	20	0.018	2	613.00

Flow Cells, Sub-micro, round aperture

Type 585.2 is designed to optimize a small sample volume by using a cylindrical sample chamber with a small diameter. The diameter of the aperture on the 585.2 series of flow cells is 2 mm. The choice of the correct 'Z' dimension for your instrument is critical as you must ensure that the light beam of the instrument passes through the sample chamber of the cell. The cell is designed with a threaded jacket which will accept M6 threaded connectors. Available 'Z' Dimensions are 8.5, 15 (cell height 35mm) and 20mm (cell height 40mm).




Catalog Number	Path Length mm		terior, n Length		Sample Chan Diameter Ø	/	'Z' Dim mm	Nominal Vol. ml	Polished Windows	Price per cell
Special Optical	Glass windov	vs, Useal	ble range	e: 320 to 2	2500 nm				M6 Threa	ded Jacket
585.2-SOG-10/Z	8.5 10	12.5	12.5	35	2	10	8.5	0.032	2	\$ 417.00
585.2-SOG-10/Z	15 10	12.5	12.5	35	2	10	15	0.032	2	417.00
585.2-SOG-10/Z	20 10	12.5	12.5	40	2	10	20	0.032	2	456.00
Spectrosil [®] Far	• UV Quartz w	vindows,	Useable	range: 1	70 to 2700 nm				M6 Thread	ed Jacket
585.2-Q-10/Z8.5	10	12.5	12.5	35	2	10	8.5	0.032	2	\$ 573.00
585.2-Q-10/Z15	10	12.5	12.5	35	2	10	15	0.032	2	573.00
585.2-O-10/Z20	10	12.5	12.5	40	2	10	20	0.032	2	613.00



Flow Cells, Round Aperture

Type 585.3 is designed to optimize a small sample volume by using a cylindrical sample chamber with a small diameter. The diameter of the aperture on the 585.3 series of flow cells is 3mm. The choice of the correct 'Z' dimension for your instrument is critical as you must ensure that the light beam of the instrument passes through the sample chamber of the cell. The cell is designed with a threaded jacket which will accept M6 threaded connectors. Available 'Z' Dimensions are 8.5, 15 (cell height 35mm) and 20mm (cell height 40mm).


'Z' Dimension	per Instrument
Manufacturer:	'Z' Dimension:
Agilent®	varies by instrument
Beckman [®]	8.5 mm
Bio-Rad [®]	8.5 mm
Eppendorf [®]	8.5 mm
GBC ^(R)	15 mm
Hewlett Packard®	15 mm
Hitachi [®]	varies by instrument
Jasco [®]	12 mm
Ocean Optics [®]	15 mm
Perkin-Elmer [®]	15 mm
Pharmacia [®]	15 mm
Shimadzu®	15 mm
StellarNet [®]	15 mm
Thermo Spectronic [®]	8.5 and 15 mm
Turner	8.5 mm
Varian [®]	20 mm

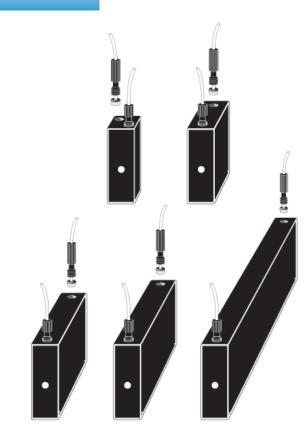


Note: Cells pictured have a Z dimension of 15mm

Catalog Number	Path Length mm		terior, n Length		Sample Cha Diameter ('Z' Dim mm	Nominal Vol. ml	Polished Windows	Price per cell
Special Optica	ıl Glass window	vs, Useal	ole range	: 320 to 2	2500 nm				M6 Thread	ded Jacket
585.3-SOG-10/2		12.5	12.5	35	3	10	15	0.070	2	\$ 403.00
585.3-SOG-20/2	Z15 20	12.5	22.5	35	3	20	15	0.140	2	621.00
585.3-SOG-40/2	Z15 40	12.5	42.5	35	3	40	15	0.280	2	914.00
585.3-SOG-50/2	Z15 50	12.5	52.5	35	3	50	15	0.350	2	1,104.00
585.3-SOG-100	/ Z15 100	12.5	102.5	35	3	100	15	0.700	2	1,710.00
(Z=8.5mm and	Z=20mm also av	ailable)								
Spectrosil [®] Fo	ar UV Quartz w	vindows,	Useable	range: 12	70 to 2700 nm				M6 Thread	ed Jacket
585.3-Q-10/Z8.	5 10	12.5	12.5	35	3	10	8.5	0.070	2	\$ 490.00
585.3-Q-20/Z8.	5 20	12.5	22.5	35	3	20	8.5	0.140	2	678.00
585.3-Q-40/Z8.	5 40	12.5	42.5	35	3	40	8.5	0.280	2	1,096.00
585.3-Q-50/Z8.	5 50	12.5	52.5	35	3	50	8.5	0.350	2	1,261.00
585.3-Q-100/Z8	3.5 100	12.5	102.5	35	3	100	8.5	0.700	2	1,750.00
	10	10.5	10.5	25	2	10	15	0.070	2	¢ 400.00
585.3-Q-10/Z15		12.5	12.5	35	3	10	15		2	\$ 490.00
585.3-Q-20/Z15		12.5	22.5	35 35	3	20	15	0.140	2	678.00
585.3-Q-40/Z15		12.5 12.5	42.5 52.5	35 35	3	40 50	15 15	0.280 0.350	2	1,096.00 1,261.00
585.3-Q-50/Z15		12.5	52.5 102.5	35 35	3	50 100	15 15	0.350	2 2	1,261.00
585.3-Q-100/Z1	15 100	12.3	102.3	33	5	100	13	0.700	Z	1,750.00
585.3-Q-10/Z20	10	12.5	12.5	40	3	10	20	0.070	2	\$ 529.00
585.3-Q-20/Z20	20	12.5	22.5	40	3	20	20	0.140	2	718.00
585.3-Q-40/Z20		12.5	42.5	40	3	40	20	0.280	2	1,130.00
585.3-Q-50/Z20	50	12.5	52.5	40	3	50	20	0.350	2	1,301.00
585.3-Q-100/Z2		12.5	102.5	40	3	100	20	0.700	2	1,797.00

The Type 31 range has path lengths 0.01mm to 5mm, with a sample compartment aperture of 15mm. Please note when filling cylindrical cells, to avoid damage to the windows, do not fill cell to the top. Leave an air gap otherwise the windows will be broken by the hydraulic pressure of the stopper being inserted. A PTFE stopper is supplied with the cell.

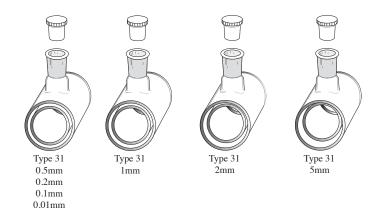

INTERNAL CROSS SECTION

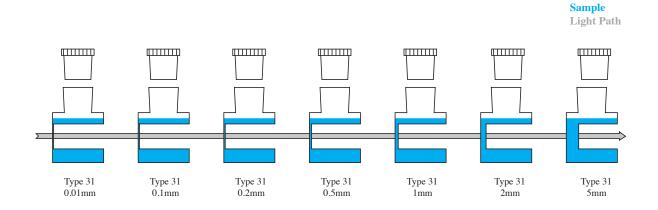

Catalog Number	Path Length mm	Exterio Diameter	/	Interio Width	or, mm Length	Nominal Vol. ml	Number of Stoppers	Polished Windows	Price per cell
Spectrosil [®]	⁹ Far UV Quartz		0	e: 170 to 270	0				L
31-Q-0.01	0.01	22	22.5	15	0.01	2.150	1	2	\$ 652.00
31-Q-0.1	0.10	22	22.5	15	0.1	2.150	1	2	443.00
31-Q-0.2	0.20	22	22.5	15	0.2	2.180	1	2	443.00
31-Q-0.5	0.50	22	22.5	15	0.5	2.220	1	2	443.00
31-Q-1	1	22	22.5	15	1	2.310	1	2	430.00
31-Q-2	2	22	22.5	15	2	2.490	1	2	425.00
31-Q-5	5	22	22.5	15	5	3.020	1	2	364.00
Infrasil [®] o	r equivalent Nea	r Infrared Q	uartz winde	ows, Useable	e range: 22	0 to 3800 nm			
31-I-0.01	0.01	22	22.5	15	0.01	2.150	1	2	\$ 678.00
31-I-0.1	0.10	22	22.5	15	0.1	2.150	1	2	464.00
31-I-0.2	0.20	22	22.5	15	0.2	2.180	1	2	464.00
31-I-0.5	0.50	22	22.5	15	0.5	2.220	1	2	464.00
31-I-1	1	22	22.5	15	1	2.310	1	2	464.00
31-I-2	2	22	22.5	15	2	2.490	1	2	430.00
31-I-5	5	22	22.5	15	5	3.020	1	2	378.00


22

The Type 31B range has path lengths 0.01mm to 5mm, with a sample compartment aperture of 13mm. The structure of this cylindrical cell allows for lower volume usage as compared to the requirements of the Type 31 cells.


Please note when filling cylindrical cells, to avoid damage to the windows, do not fill cell to the top. Leave an air gap otherwise the windows will be broken by the hydraulic pressure of the stopper being inserted. PTFE stoppers are supplied with these cells as shown in drawing below.

0.2mm 0.5mm

The tubulations are 7mm OD x 5mm ID.

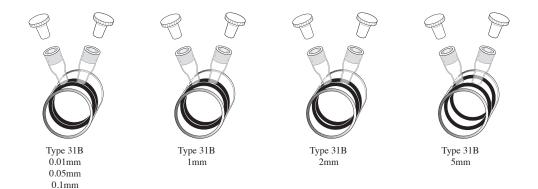

INTERNAL CROSS SECTION

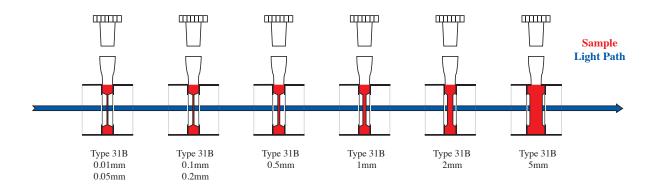

Catalog Number	Path Length mm	Exterio Diameter	or, mm Length	Interio Width	or, mm Length	Nominal Vol. ml	Number of Stoppers	Polished Windows	Price per cell
Spectrosil [®]	Far UV Quartz	windows, Us	seable range	e: 170 to 270	00 nm				
31B-Q-0.01	0.01	22	22.5	13	0.01	0.140	2	2	\$ 907.00
31B-Q-0.05	0.05	22	22.5	13	0.05	0.151	2	2	786.00
31B-Q-0.1	0.1	22	22.5	13	0.1	0.165	2	2	523.00
31B-Q-0.2	0.2	22	22.5	13	0.2	0.194	2	2	523.00
31B-Q-0.5	0.5	22	22.5	13	0.5	0.278	2	2	523.00
31B-Q-1	1	22	22.5	13	1	0.420	2	2	523.00
31B-Q-2	2	22	22.5	13	2	0.703	2	2	496.00
31B-Q-5	5	22	22.5	13	5	1.552	2	2	496.00

Type 32 and 34 cells range from path lengths of 10mm to 100mm with an internal diameter of 19mm. Other diameters available from 5mm ID to 30mm ID on request. Please note when filling cylindrical cells, to avoid damage to the windows, do not fill cell to the top. Leave an air gap otherwise the windows will be broken by the hydraulic pressure of the stopper being inserted. PTFE stoppers are supplied with each cell as pictured.

Adaptor for IR instruments


These adaptors will allow the use of 10 through 100mm cylindrical cells with Infrared instruments. All come complete with a 3" x 2" back plate to fit standard Infrared instrument holders. The adaptors are constructed of black anodized aluminum.

Cat. No	Description	Price, each	
CH-34/100	Cylindrical 50 to	100mm \$242.00	



Catalog Number	Path Length mm	Exterio Diameter	/	Interio Width	or, mm Length	Nominal Vol. ml	Number of Stoppers	Polished Windows	Price per cell
Special Opt	ical Glass windo	ows, Useable	range: 320	to 2500 nm					
32-SOG-10	10	22	12.5	19	10	2.800	1	2	\$ 124.00
32-SOG-20	20	22	22.5	19	20	5.600	1	2	137.00
34-SOG-50	50	22	52.5	19	50	14.100	2	2	150.00
34-SOG-100	100	22	102.5	19	100	28.200	2	2	164.00
Spectrosil [®]	Far UV Quartz	windows, Us	seable rang	e: 170 to 270	00 nm				
32-Q-10	10	22	12.5	19	10	2.800	1	2	\$ 169.00
32-Q-20	20	22	22.5	19	20	5.600	1	2	181.00
34-Q-50	50	22	52.5	19	50	14.100	2	2	242.00
34-Q-100	100	22	102.5	19	100	28.200	2	2	268.00
Infrasil [®] or	· equivalent Nea	r Infrared Q	uartz wind	ows, Useabl	e range: 22	20 to 3800 nm			
32-I-10	10	22	12.5	19	10	2.800	1	2	\$ 189.00
32-I-20	20	22	22.5	19	20	5.600	1	2	203.00
34-I-50	50	22	52.5	19	50	14.100	2	2	260.00
34-I-100	100	22	102.5	19	100	28.200	2	2	294.00




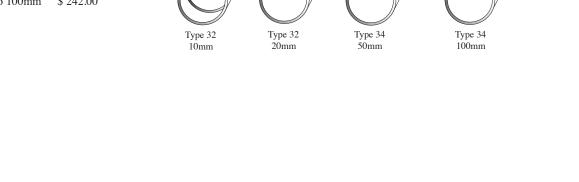

Type 35 cells have a large aperture for use where a large surface irradiation area is required in either short or long path lengths. A PTFE stopper is supplied with the type 35 cells.

Note: Liquids do not compress like gases and therefore excessive pressure applied to the stopper of an overfilled cell will cause damage to the cell.

Although a tube can be attached to virtually any cell, the Type 37 design already incorporates a tube so that the cell can be attached to other apparatus, or sealed off with a flame, depending on the application. The tube can also be specified as a quartz to pyrex graded seal. Where this is required, call us for a quotation.

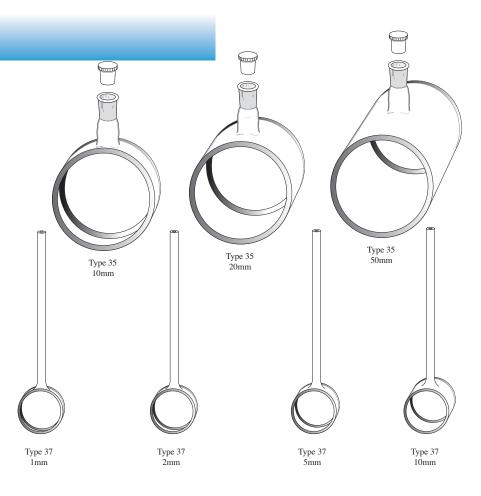


Catalog Number	Path Length mm	Exterior Diameter	·	Interio Diameter	· ·	Nominal Vol. ml	Number of Stoppers	Polished Windows	Price per cell
Optical Pyr	rex windows, Use	able range: 3	320 to 2500	nm					
35-PX-10	10	50	12.5	47	10	17.00	1	2	\$ 228.00
35-PX-20	20	50	22.5	47	20	35.00	1	2	242.00
35-PX-50	50	50	52.5	47	50	86.00	1	2	273.00
Spectrosil [®]	Far UV Quartz	windows, Us	eable range	e: 170 to 270	00 nm				
35-Q-10	10	50	12.5	47	10	17.00	1	2	\$ 438.00
35-Q-20	20	50	22.5	47	20	35.00	1	2	529.00
35-Q-50	50	50	52.5	47	50	86.00	1	2	608.00
Optical Pyr	rex windows, Use	eable range: 3	320 to 2500	nm			Tube Length mm		
37-PX-1	1	22	3.5	19	1	0.28	70	2	\$ 124.00
37-PX-2	2	22	4.5	19	2	0.56	70	2	124.00
37-PX-5	5	22	7.5	19	5	1.40	70	2	124.00
37-PX-10	10	22	12.5	19	10	2.80	70	2	124.00
Spectrosil [®]	^y Far UV Quartz	windows, Us	eable range	e: 170 to 270	00 nm				
37-Q-1	1	22	3.5	19	1	0.28	70	2	\$ 233.00
37-Q-2	2	22	4.5	19	2	0.56	70	2	233.00
37-Q-5	5	22	7.5	19	5	1.40	70	2	221.00
37-Q-10	10	22	12.5	19	10	2.80	70	2	216.00
Infrasil [®] o	or equivalent Nea	r Infrared Q	uartz windo	ows, Useable	e range: 22	0 to 3800 nm			
37-I-1	1	22	3.5	19	1	0.28	70	2	\$ 247.00
37-I-2	2	22	4.5	19	2	0.56	70	2	247.00
37-I-5	5	22	7.5	19	5	1.40	70	2	242.00
37-I-10	10	22	12.5	19	10	2.80	70	2	228.00

Fluorometer Cells, Standard, Semi-Micro, Micro

Standard rectangular fluorometer cells Type 3 are the most popular cells for general fluorescence measurement and they will fit virtually every fluorometer instrument. All four windows and the bottom of the cell are polished and because of the technique used in their construction they remain extremely flat, right to the very edge of the window, allowing them to be used in some laser applications. Each cell is supplied with a PTFE cover.

Smaller volume Types 9F, 17F & 18F are identical to the Type 9 Semi-Micro, Type 17 Short-Micro and Type 18 Micro range, except that all cell walls and the base are optically polished.

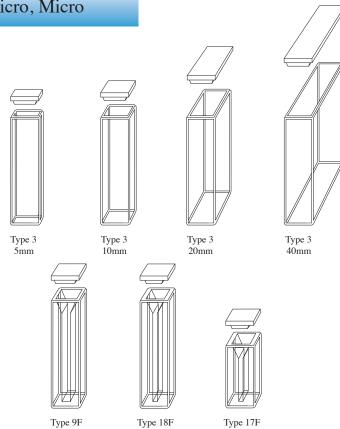

10mm

10mm

Catalog	Path Length	Ex	xterior, m	m	Interio	or, mm	Nominal	Polish	ed	Price
Number	mm	Width	Length	Height	Width	Length	Vol. ml	Windows	Base	per cell
Optical Glass	s windows, Useable	e range: 3	34 to 250	0 nm						
3-G-5	5	12.5	7.5	45	10	5	1.700	4	1	\$ 103.00
3-G-10	10	12.5	12.5	45	10	10	3.500	4	1	94.00
3-G-20	20	12.5	22.5	45	10	20	7.000	4	1	164.00
3-G-40	40	12.5	42.5	45	10	40	14.00	4	1	208.00
Special Optic	al Glass windows,	Useable 1	range: 320) to 2500 1	ım					
3-SOG-10	10	12.5	12.5	45	10	10	3.500	4	1	\$ 137.00
Spectrosil [®] I	Far UV Quartz win	dows, Us	eable rang	ge: 170 to	2700 nm					
3-Q-5	5	12.5	7.5	45	10	5	1.700	4	1	\$ 255.00
3-Q-10	10	12.5	12.5	45	10	10	3.500	4	1	181.00
3-Q-20	20	12.5	22.5	45	10	20	7.000	4	1	299.00
3-Q-40	40	12.5	42.5	45	10	40	14.00	4	1	351.00
Infrasil [®] or	equivalent Near In	frared Qu	uartz wind	lows, Use	able range: 2	220 to 3800 n	m			
3-I-5	5	12.5	7.5	45	10	5	1.700	4	1	\$ 260.00
3-I-10	10	12.5	12.5	45	10	10	3.500	4	1	242.00
3-I-20	20	12.5	22.5	45	10	20	7.000	4	1	346.00
3-I-40	40	12.5	42.5	45	10	40	14.00	4	1	403.00
Spectrosil [®] I	Far UV Quartz win	dows, Us	eable rang	ge: 170 to	2700 nm			Semi-l	Micro,	Ope8Top
9F-Q-10	10	12.5	12.5	45	4	10	1.400	4	1	\$ 320.00
								M	icro, (Open Top
17F-Q-10	10	12.5	12.5	25	2	10	0.400	4	1	\$ 281.00
18F-Q-10	10	12.5	12.5	45	2	10	0.700	4	1	351.00


10mm

26

Fluorometer Cells, Standard, Semi-Micro, Micro with stopper

Type 23 cells are identical to Standard Rectangular fluorometer cells as listed on page 26 except that instead of having an open rectangular hole at the top, a block is fused to the top of the cell with a ground hole to receive a PTFE stopper. This provides a more suitable seal for volatile liquids. A PTFE stopper is supplied with each cell. Two stoppers are supplied with cells which have a path length of 40mm or greater.

Smaller volume Types 29F and 28F are identical to Types 9F Semi-Micro and Type 18F Micro cell ranges except they have a stopper instead of a lid. A PTFE stopper is supplied with each cell.

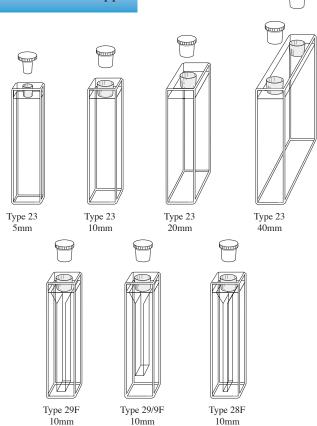

-

Catalog	Path Length	Ex	terior, m	m	Interi	or, mm	Nominal	Polish	ed	Price
Number	mm	Width	Length	Height	Width	Length	Vol. ml	Windows	Base	per cell
Optical Glass	s windows, Useabl	e range: 3	34 to 250	00 nm						
23-G-5	5	12.5	7.5	48	10	5	1.700	4	1	\$ 150.00
23-G-10	10	12.5	12.5	48	10	10	3.500	4	1	129.00
23-G-20	20	12.5	22.5	48	10	20	7.000	4	1	242.00
23-G-40	40	12.5	42.5	48	10	40	14.00	4	1	286.00
Special Optic	cal Glass windows,	, Useable 1	range: 32	0 to 2500 i	nm					
23-SOG-10	10	12.5	12.5	48	10	10	3.500	4	1	\$ 189.00
Spectrosil [®] 1	Far UV Quartz wit	ndows, Us	eable ran	ge: 170 to	2700 nm					
23-Q-5	5	12.5	7.5	48	10	5	1.700	4	1	\$ 299.00
23-Q-10	10	12.5	12.5	48	10	10	3.500	4	1	242.00
23-Q-20	20	12.5	22.5	48	10	20	7.000	4	1	378.00
23-Q-40	40	12.5	42.5	48	10	40	14.00	4	1	464.00
Infrasil [®] or	equivalent Near Ir	nfrared Qi	iartz win	dows, Use	able range:	220 to 3800 n	ım			
23-I-5	5	12.5	7.5	48	10	5	1.700	4	1	\$ 346.00
23-I-10	10	12.5	12.5	48	10	10	3.500	4	1	325.00
23-I-20	20	12.5	22.5	48	10	20	7.000	4	1	417.00
23-I-40	40	12.5	42.5	48	10	40	14.00	4	1	495.00
Spectrosil [®] 1	Far UV Quartz wit	ndows, Us	eable ran	ge: 170 to	2700 nm			Sem	i-Micr	o, Stopprr
29F-Q-10	10	12.5	12.5	48	4	10	1.400	4	1	\$ 386.00
29/9F-Q-10	10	12.5	12.5	48	4	10	1.160	4	1	430.00
									Micro	o, Stopper
28F-Q-10	10	12.5	12.5	48	2	10	0.700	4	1	\$ 403.00




Fluorometer Micro Square Cells, Adaptors

This range of Micro Fluorometer cells has been specifically designed by Starna for use with the FCA range of adaptors listed below. The adaptor locates the square cross section cell accurately in the center of the instrument cell compartment. Therefore, for a given path length, the cell is optimized for minimum volume and maximum sensitivity. All four windows and the base of the cells are polished and are available with either an open top or with a PTFE stopper supplied with the type 23- cells.

The FCA adaptors have exterior dimensions of a standard fluorometer cell (12.5mm x 12.5mm) and interior dimensions to accommodate the square micro fluorometer cells listed. Apertures are machined into the adaptor wall to allow for transmission of the excitation and emission energy for instrument 'Z' dimensions of 8.5, 15 or 20mm. The aluminum adaptors are black anodized to eliminate stray light. The correct adaptor for a particular cell can be identified from the listing below.

Adaptor	Adaptor Aperture Dimensions										
	Width	Height	Number of								
	mm	mm	Openings								
FCA3.30	2.5	11	4								
FCA3	2.5	21	4								
FCA4	3.5	21	4								
FCA5	4.5	21	4								
	From the bottom of the adaptor to the bottom of the opening is 5.8 mm for all adaptors.										

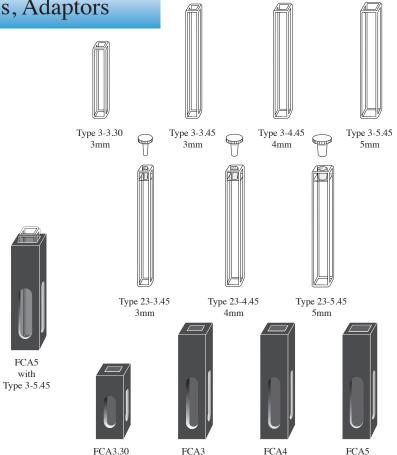

Catalog Number	Path Length mm		xterior, n Length		Interio Width	or, mm Length	Adapter	Nominal Vol. ml	Polished Windows	Price per cell
Special Optical Glass windows, Useable range: 320 to 2500 nm Open Top										
3-3.30-SOG-3	3	5.5	5.5	30	3	3	FCA3.30	0.225	5	\$ 137.00
3-3.45-SOG-3	3	5.5	5.5	45	3	3	FCA3	0.315	5	137.00
3-4.45-SOG-4	4	6.5	6.5	45	4	4	FCA4	0.560	5	137.00
3-5.45-SOG-5	5	7.5	7.5	45	5	5	FCA5	0.875	5	137.00
										Stopper
23-3.45-SOG-3	3	5.5	5.5	48	3	3	FCA3	0.315	5	\$ 164.00
23-4.45-SOG-4	4	6.5	6.5	48	4	4	FCA4	0.560	5	164.00
23-5.45-SOG-5	5 5	7.5	7.5	48	5	5	FCA5	0.875	5	164.00
Spectrosil [®] F	ar UV Quartz	windows,	Useable 1	range: 170	to 2700 nm	!				Open Top
3-3.30-Q-3	3	5.5	5.5	30	3	3	FCA3.30	0.225	5	\$ 307.00
3-3.45-Q-3	3	5.5	5.5	45	3	3	FCA3	0.315	5	307.00
3-4.45-Q-4	4	6.5	6.5	45	4	4	FCA4	0.560	5	307.00
3-5.45-Q-5	5	7.5	7.5	45	5	5	FCA5	0.875	5	307.00
										Stopper
23-3.45-Q-3	3	5.5	5.5	48	3	3	FCA3	0.315	5	\$ 320.00
23-4.45-Q-4	4	6.5	6.5	48	4	4	FCA4	0.560	5	320.00
23-5.45-Q-5	5	7.5	7.5	48	5	5	FCA5	0.875	5	320.00

Adapters for the above series of Micro cells

Catalog		Ex	xterior, m	m	Interi	or, mm		Price
Number	Description	Width	Length	Height	Width	Length	Fit cell types	per adapter
FCA3.30	Adapter	12.5	12.5	30	5.55	5.55	3-3.30	\$ 155.00
FCA3	Adapter	12.5	12.5	45	5.55	5.55	3-3, 23-3	155.00
FCA4	Adapter	12.5	12.5	45	6.55	6.55	3-4, 23-4	155.00
FCA5	Adapter	12.5	12.5	45	7.55	7.55	3-5, 23-5	155.00



Fluorometer Sub-Micro

Sub-Micro cells, Type 16F, while retaining the exterior dimensions of a standard cell, are designed for the measurement of very small samples with volumes ranging from 10μ l to 160μ l. The fluorometer Sub-Micro cells have a third window at 90° to the incident light to allow emission energy to be measured. Special adaptors are not required for this cell.

The entrance to the sample compartment is hemispherical, being designed without sharp corners to eliminate potential loss of sample by capillary action. The amount of sample required to fill the sample chamber is reduced to an absolute minimum, typically fifteen to twenty per cent greater than the absolute sample chamber volume. The sample is easily inserted into and retrieved from the cell by careful use of a pipette or syringe.

The sample compartment in the Sub-micro range of cells optimizes the alignment of the sample for maximum sensitivity in the fluorometer. The correct 'Z' dimension needs to be selected for this range, typically 8.5 or 15 mm. Open top cells are supplied with a tight sealing polypropylene cap, as well as a PTFE lid.


'Z' Dimensions for some fluorometers

Manufacturer:	'Z' Dimension:
	varies by instrument
Jasco [®]	18mm
Molecular Devices	® 15mm
Ocean Optics [®]	15mm
Perkin-Elmer [®]	15mm
Pharmacia [®]	15mm
PTI (Photon Techn	ology) [®] 15mm
Shimadzu [®]	15mm
SLM/Spectronics®) 15mm
Spectra Max [®]	15mm
Spex®	15mm
StellarNet®	15mm
TSS [®]	15mm
Varian [®]	20mm







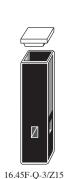








16.160F-Q-10/Z15

16.10F-Q-10/Z15

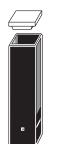

16.10F-Q-10/Z8.5



И

16.12F-Q-1.5/Z15

16.40F-Q-10/Z8.5



16.45F-Q-3/Z8.5

16.50F-Q-10/Z8.5

16.100F-Q-10/Z8.5

16.160F-Q-10/Z8.5

Catalog Number	Path Length mm		terior, n Length			e Chamb Length	· · · · · · · · · · · · · · · · · · ·	'Z' Dim mm	Nominal Vol. ml	Polished Windows	Price per cell
Spectrosil [®] Far	UV Quartz v	vindows,	Useable	range: 1	70 to 270	00 nm					Open Top
16.10F-Q-10/Z15	10	12.5	12.5	45	1	10	1	15	0.010	3	\$ 568.00
16.10F-Q-10/Z8.5	10	12.5	12.5	45	1	10	1	8.5	0.010	3	568.00
16.12F-Q-1.5/Z15	1.5	12.5	12.5	45	1.5	1.5	5	15	0.012	3	568.00
16.12F-Q-1.5/Z8.4	5 1.5	12.5	12.5	45	1.5	1.5	5	8.5	0.012	3	568.00
16.40F-Q-10/Z15	10	12.5	12.5	45	2	10	2	15	0.040	3	542.00
16.40F-Q-10/Z8.5	10	12.5	12.5	45	2	10	2	8.5	0.040	3	542.00
16.45F-Q-3/Z15	3	12.5	12.5	45	3	3	5	15	0.045	3	529.00
16.45F-Q-3/Z8.5	3	12.5	12.5	45	3	3	5	8.5	0.045	3	529.00
16.50F-Q-10/Z15	10	12.5	12.5	45	2	10	2.5	15	0.050	3	516.00
16.50F-Q-10/Z8.5	10	12.5	12.5	45	2	10	2.5	8.5	0.050	3	516.00
16.100F-Q-10/Z1	5 10	12.5	12.5	45	2	10	5	15	0.100	3	503.00
16.100F-Q-10/Z8.	5 10	12.5	12.5	45	2	10	5	8.5	0.100	3	503.00
16.160F-Q-10/Z1	5 10	12.5	12.5	45	2	10	8	15	0.160	3	503.00
16.160F-Q-10/Z8.	5 10	12.5	12.5	45	2	10	8	8.5	0.160	3	503.00

www.starnacells.com




16.100F-Q-10/Z15





Fluorometer, Triangular, Dual Path Length

Triangular fluorometer cells Type 4, with open top or Type 24 with PTFE stopper, are useful for front surface fluorescence measurement. This technique is employed when the sample is too dense for the transmission of excitation or emission energy. Analysis is made by exciting the face of the sample at an angle of 45° to both the excitation source and the emission detector.

Dual path length cells Types 52 and 53 allow for both fluorescence and absorbance measurements to be made with small volumes with a choice of two path lengths, dependent on the orientation of the cell in the cell holder. All cells have a 10mm path length in one direction and a smaller path length at 90°. All four windows and the base are clear. Type 53 is supplied with a PTFE stopper.

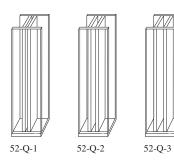

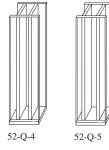

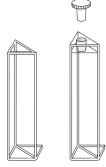

4-SB-Q-10 24-SB-Q-10

					53-Q-1	53-Q-2	53-Q-3	53-Q-4	53-Q-5
Catalog Number	Path Length mm	E: Width	xterior, n Length		Interi Width	ior, mm Length	Nominal Vol. ml	Polished Windows	Price per cell
Spectrosil [®] F	ar UV Quartz wit	ndows, Us	eable ran	ge: 170 to	o 2700 nm			Triangular,	Open Top
4-SB-Q-10	Surface	12.5	12.5	45	10	10	1.700	3	\$ 299.00
								Triangulo	ır, Stopper
24-SB-Q-10	Surface	12.5	12.5	48	10	10	1.700	3	\$ 378.00


M

Catalog Number	Path Length mm	E Width	xterior, n Length	nm Height	Interi Width	or, mm Length	Nominal Vol. ml	Polished Windows	Price per cell
Spectrosil [®]	Far UV Quartz wit	ndows.Us	eable ran	ge: 170 to	2700 nm	8		Dual Path,	Open Top
52-0-1	1 and 10	12.5	12.5	45	1	10	0.400	5	\$ 281.00
52-Q-2	2 and 10	12.5	12.5	45	2	10	0.800	5	281.00
52-Q-3	3 and 10	12.5	12.5	45	3	10	1.200	5	281.00
52-Q-4	4 and 10	12.5	12.5	45	4	10	1.600	5	281.00
52-Q-5	5 and 10	12.5	12.5	45	5	10	2.000	5	281.00
-								Dual Path, wi	th Stopper
53-Q-1	1 and 10	12.5	12.5	48	1	10	0.400	5	\$ 359.00
53-0-2	2 and 10	12.5	12.5	48	2	10	0.800	5	359.00
53-Q-3	3 and 10	12.5	12.5	48	3	10	1.200	5	359.00
53-0-4	4 and 10	12.5	12.5	48	4	10	1.600	5	359.00
53-Q-5	5 and 10	12.5	12.5	48	5	10	2.000	5	359.00



Fluorometer, Tandem or Divided Rectangular


Types 55, 56, 57 and 58 are all cells with two chambers which can be used either in parallel or series and with the exception of Type 55-1, have all walls polished. Types 57 and 58 have the facility for the intentional mixing of two liquids within the cell itself. PTFE stoppers supplied for type 56 and type 58 as pictured.

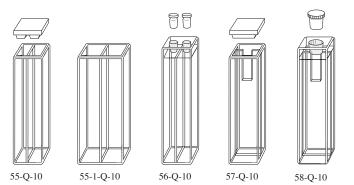


Catalog Number	Path Length mm	E Width	xterior, n Length		Inte Width	rior, mm Length	Nominal Vol. ml	Polished Windows	Price per cell
Spectrosil [®] Far UV Quartz windows, Useable range: 170 to 2700 nm Tandem or Divided Rectangular, Open top									
55-Q-10	2 x 10	12.5	12.5	45	10	2 x 4.375	2 x 1.500	6	\$ 281.00
55-1-Q-10	2 x 10	12.5	23.75	45	10	2 x10	2 x 3.500	3	325.00
						Тат	dem or Divided Re	ectangular, wit	h Stoppers
56-Q-10	2 x 10	12.5	12.5	48	10	2 x 4.375	2 x 1.500	6	\$ 430.00
						Mixing	, Tandem or Divide	ed Rectangular	, Open top
57-Q-10	2 x 10	12.5	12.5	45	10	2 x 4.375	2 x 1.000	6	\$ 364.00
						Mixing, Tar	dem or Divided Re	ectangular, wit	h Stoppers
58-Q-10	2 x 10	12.5	12.5	48	10	2 x 4.375	2 x 1.000	6	\$ 417.00

Fluorometer Screw Cap Cells, Closed and Septum

These fluorometer cells with the suffix GL14 are able to be used under anaerobic conditions. The cells with the GL14 suffix have a standard GL14 threaded top and cap which can be either a septum cap, GL14-S, or a plain closed cap, GL14-C.

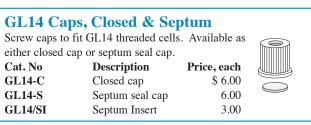

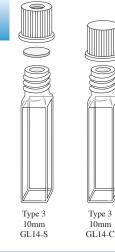

Type 3

Type 9F

10mm

GL14-S



GL14-C

Catalog Number	Patl Description	h Length mm		xterior, r Length			or, mm Length	Nominal Vol. ml	Polisho Windows		Price se each
	lass windows, Useable rang	e: 334 t		U	8		8				
· ·	Closed Cap, Fluorometer	10	12.5	12.5	68	10	10	3.500	4	1	\$ 273.00
Spectrosil [®] Far UV Quartz windows, Useable range: 170 to 2700 nm											
3-Q-10-GL14-C	Closed Cap, Fluorometer	10	12.5	12.5	68	10	10	3.500	4	1	\$ 359.00
9F-Q-10-GL14-C	Closed Cap, Semi-Micro,F	10	12.5	12.5	68	4	10	1.400	4	1	469.00
Special Optical G	lass windows, Useable rang	ge: 334 t	o 2500 i	nm							
3-SOG-10-GL14-S	Septum Cap, Fluorometer	10	12.5	12.5	68	10	10	3.500	4	1	\$ 273.00
Spectrosil [®] Far UV Quartz windows, Useable range: 170 to 2700 nm											
3-Q-10-GL14-S	Septum Cap, Fluorometer	10	12.5	12.5	68	10	10	3.500	4	1	\$ 359.00
9F-Q-10-GL14-S	Septum Cap, Semi-Micro,F	10	12.5	12.5	68	4	10	1.400	4	1	469.00

Fluorometer Flow Through

Type 583-F range of fluorometer flow cells is designed to fit most standard fluorometers. The inlet/outlet tubes are connected with a threaded M6 connector with PTFE tubing attached. A pair of M6 connectors is supplied with the type 583 cell. The cells have three clear windows. The excitation beam passes through the cell, which also allows the cell to be used for absorption spectroscopy as well. The emmission window is at 90 degrees to the excitation windows.

The cells are available in several configurations offering a range of sample volumes. The 'Z' dimension is the distance from the bottom of the cell holder to the center of the excitation beam. It is important to match the cell 'Z' dimension with that of your instrument.









583.2.2F-Q-10

583.4F-Q-10

Manufacturer:	'Z' Dimension:
Hitachi®	varies by instrument
Jasco [®]	18mm
	®

'Z' Dimensions for some fluorometers

Jasco	18mm
Molecular Devices®	15mm
Ocean Optics [®]	15mm
Perkin-Elmer [®]	15mm
Pharmacia [®]	15mm
PTI (Photon Technology) [®]	15mm
Shimadzu®	15mm
SLM/Spectronics [®]	15mm
Spectra Max [®]	15mm
Spex ^(R)	15mm
StellarNet [®]	15mm
TSS [®]	15mm
Varian [®]	20mm



Flow cell Tubulation dimensions:								
Туре	OD	ID	Length					
Non-Threaded:	3.5mm	2.2mm	10mm					
M6	1.6mm	1mm	400mm					

Extra M6 tubing sets may be ordered.

Each set has one inlet and one outlet tube each with an M6 connector. The tubing supplied is approximately 1.5 meters long.

Catalog#	Description	Price/set
M6-Set	M6 Inlet/Outlet tubes	\$ 111.00


Catalog Number	Path Length mm		xterior, 1 Length		1		ber, mm Height	'Z' Dim mm	Nominal Vol ml	Polished Windows	Price per cell
Spectrosil [®] Fa	r UV Quartz	windows	s, Useabl	e range: 1	170 to 270	0 nm					
583.2.2F-Q-10/Z	15 10	12.5	12.5	35	2	10	2	15	0.040	3	\$ 908.00
583.2.2F-Q-10/Z	2 20 10	12.5	12.5	35	2	10	2	20	0.040	3	908.00
583.4.2F-Q-10/Z	2 15 10	12.5	12.5	35	4	10	2	15	0.080	3	940.00
583.4.2F-Q-10/Z	2 20 10	12.5	12.5	35	4	10	2	20	0.080	3	940.00
583.1F-Q-1/Z15	1	12.5	12.5	35	1	1	11	15	0.011	3	940.00
583.3F-Q-3/Z15	3	12.5	12.5	35	3	3	11	15	0.100	3	953.00
583.4F-Q-10/Z1	5 10	12.5	12.5	35	4	10	11	15	0.440	3	992.00




Fluorometer Flow Through

Fluorometer cells Type 46F and 47F (Semi-Micro) are for inline continuous flow type applications, usually measuring at one wavelength. Samples normally flow from the bottom to the top.

Types 71F and 72F (Semi-Micro) fluorescent cells have long windows and will fit normal cell holders. They can be used with instruments of any 'Z' dimension and are particularly suitable for any application where sample volume is not critical.

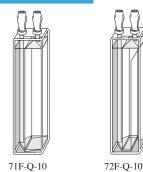

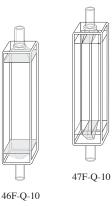

Catalog Number	Path Length mm		xterior, 1 Length	nm Height	-	e Chamb Length	oer, mm Height	'Z' Dim mm	Nominal Vol ml	Polished Windows	Price per cell
Spectrosil [®] F	ar UV Quartz v	windows	, Useabl	e range:	170 to 270	0 nm				Standar	d in-Line
46F-Q-10	10	12.5	12.5	45	10	10	40	all	4.000	4 Semi-Micr	\$ 660.00 o in-Line
47F-Q-10	10	12.5	12.5	45	4	10	40	all	1.600 Rect	4 tangular lon	\$ 769.00 g window
71F-Q-10	10	12.5	12.5	48	7	10	40	all Recta	3.000 ngular long	3 g window Se	\$ 686.00 <i>mi-Micro</i>
72F-Q-10	10	12.5	12.5	48	4	10	40	all	1.800	3	\$ 712.00

Fluoresence Reference Set

Molecular fluorescence spectroscopy is a sensitive and often selective technique. Unlike absorption spectrophotometry, it is not an absolute technique: instruments therefore require calibration before every series of measurements. This may be achieved using a stable reference material, which should absorb and emit at similar wavelengths to the samples of interest. Use of the general purpose fluorescent reference material set type 6BF enables the day to day stability of fluoresence instruments to be verified.

The 6BF reference materials are not standards with absolute values, but a set of six relatively stable fluorescent materials in a polymethylmethacrylate matrix with which the stability of the instrument can be monitored. Four blocks exhibit broad band spectra which cover the normally used UV and visible region of the spectrum. Two blocks contain materials suitable for illustrating the selectivity of the technique as well as checking instrument resolution and wavelength calibration.



Acrylic polymer references for Fluoresence

Stability:	no degradation, no evaporation
Safety:	no chemicals to mix
Robust:	unbreakable, easy to store and use
For more inform	ation call for our 6BF technical reference

Catalog Number 6BF


Description Fluorescence Reference Set, 6 polymer blocks Price per Set \$ 1,221.00




Dye Laser Cells


Dye Laser Cells must be made extremely accurate with a surface flatness which extends to the edge of the cell. Many of the fluorometer type of cells illustrated in this catalog may be used for laser applications. However, to ensure their optical flatness, fluorometer cells must be polished to more exacting tolerances. The Dye Laser Cells have a window polish of 20/10 Scratch/Dig and are given the prefix **DLC**.


Specifically designed cells for laser applications are listed below. The DLC-3-Q-10 is supplied with a PTFE lid while the other dye laser cells shown in the top row are supplied with a PTFE stopper. The DLC300-Q-20 and DLC300-Q-40 are supplied without a lid.

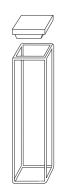




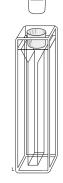

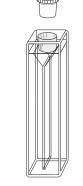

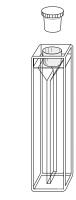

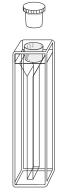

DLC-3-Q-10


DLC-23-Q-10 DLC-29F-Q-10

DLC-28.05F-Q-10

DLC-28.1F-Q-10 E

DLC-28F-Q-10




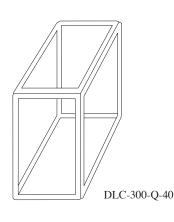



Catalog	Path Length	Ex	cterior, m	ım	Interi	or, mm	Nominal	Polish	ed	Price
Number	mm	Width	Length	Height	Width	Length	Vol. ml	Windows	Base	per cell
Spectrosil [®] Fa	r UV Quartz win	idows, Use	eable ran	ge: 170 to	2700 nm					
DLC-3-Q-10	10	12.5	12.5	45	10	10	3.500	4	1	\$ 372.00
DLC-23-Q-10	10	12.5	12.5	48	10	10	3.500	4	1	430.00
DLC-29F-Q-10	10	12.5	12.5	48	4	10	1.400	4	1	503.00
DLC-28.05F-Q-2	10 10	12.5	12.5	48	0.5	10	0.175	4	1	587.00
DLC-28.1F-Q-1	0 10	12.5	12.5	48	1	10	0.350	4	1	608.00
DLC-28F-Q-10	10	12.5	12.5	48	2	10	0.700	4	1	561.00
DLC-300-Q-20	20	26	26	40	20	20	12.000	4	no S	5 1,170.00
DLC-300-Q-40	40	26	46	40	20	40	24.000	4	no	1,397.00


Cells with Quartz to Pyrex Graded Seal

Quartz to Pyrex graded seals can be fused to most rectangular quartz cells, allowing them to be joined to other borosilicate apparatus. The price indicated is for the graded seal, the cost of fusing to a cell plus the cost of the chosen cell. The dimensions of the normal graded seals used are OD 10 mm, ID 8 mm, Length 70mm.

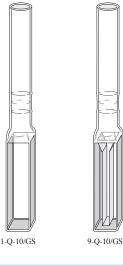

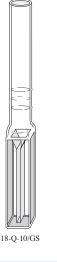

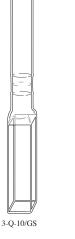


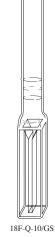

Catalog Path Length Exterior, mm Interior, mm Nominal Polished Price Number mm Width Length Height Width Length Vol. ml Windows per cell Spectrosil[®] Far UV Quartz windows, Useable range: 170 to 2700 nm 1-O-10/GS 10 12.5 12.5 115 10 3.500 2 \$ 299.00 10 9-0-10/GS 2 10 12.5 12.5 115 4 10 1.400 364.00 18-Q-10/GS 10 12.5 12.5 115 2 10 0.700 2 372.00 3-0-10/GS 10 12.5 12.5 115 10 10 3.500 5 372.00 9F-Q-10/GS 10 12.5 12.5 115 4 10 1.400 5 508.00 18F-Q-10/GS 10 12.5 2 0.700 5 12.5 115 10 542.00

Cells with Straight Bore Tube

Straight bore tubes in a range of sizes can be attached to most cells. Quartz tubes can be fused to quartz cells but for Glass materials, the tube and the cell need to be fabricated from borosilicate glass. Prices quoted are for the cost of the tube and fusing it to the cell plus the original cell cost. The dimensions of the standard tubes supplied are OD 10 mm, ID 8 mm, Length 70mm.

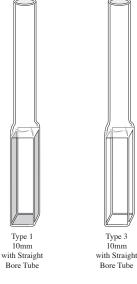


9F-Q-10/GS

Catalog Number	Path Length		terior, m Length		Interi Width	or, mm Length	Nominal Vol. ml	Polished Windows	Price per cell
Inullibel	mm	wium	Length	Height	wiuui	Length	V01. IIII	vv muows	per cen
Special Optica	l Glass windows,	Useable r	ange: 32	0 to 2500 i	nm				
1-SOG-10/SBT	10	12.5	12.5	115	10	10	3.500	2	\$ 228.00
Spectrosil [®] Fo	ır UV Quartz wir	idows, Use	eable ran	ge: 170 to	2700 nm				
1-Q-10/SBT	10	12.5	12.5	115	10	10	3.500	2	\$ 286.00
3-Q-10/SBT	10	12.5	12.5	115	10	10	3.500	5	359.00

Colorimeter Cells


Colorimeter cells are fully fused glass cells for use with colorimeters. Each cell has two polished windows and an open top with no cover.

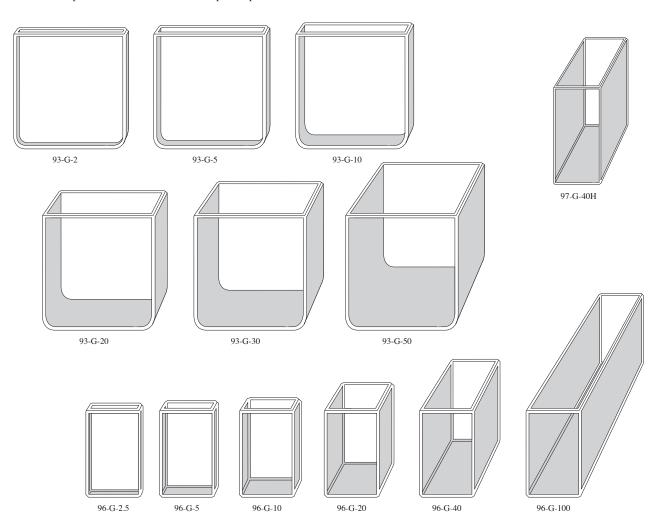



Catalog Number	Path Length mm	E: Width	xterior, m Length		Interio Width	or, mm Length	Nominal Vol. ml	Polished Windows	Price per cell
For use wit	h: Hunter Colori	imeter, AC	'S, Data C	Color and o	others				
93-G-2	2	55	2	56	51	2	25	2	\$ 260.00
93-G-5	5	55	9	56	51	5	12.5	2	209.00
93-G-10	10	55	14	56	51	10	25	2	189.00
93-G-20	20	55	24	56	51	20	50	2	233.00
93-G-30	30	55	34	56	51	30	75	2	325.00
93-G-50	50	55	54	56	51	50	125	2	378.00
For use wit	h: CP method an	d general	colorime	tric work					
97-G-40H	40	20	44	45	16	40	25	2	\$ 103.00
For use wit	h: Absorptiomete	rs							
96-G-2.5	2.5	28	8.5	40	24	2.5	1.8	2	\$ 55.00
96-G-5	5	28	11	40	24	5	3.6	2	55.00
96-G-10	10	28	16	40	24	10	7.2	2	55.00
96-G-20	20	28	26	40	24	20	14	2	77.00
96-G-40	40	28	46	40	24	40	28	2	103.00
96-G-100	100	28	106	40	24	100	70	2	164.00


36


Cells for use with Magnetic Stirrers

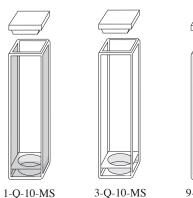


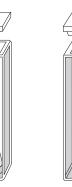

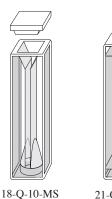

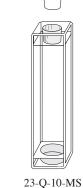

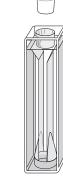

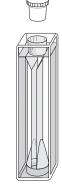

1-Q-10-MS

9-Q-10-MS

This range of spectrophotometer cells is designed specifically for use with magnetic stir bar type stirrers. The circular recess in the bottom both contains the stir bar and increases the stirring efficiency. Cells are supplied with a PTFE lid or stopper as pictured. 21-Q-10-MS

29-Q-10-MS


28-Q-10-MS

Catalog		Windows		Nominal	Price
Number	Description	Clear	Material	Volume	per cell
1-Q-10-MS	Stirring Cell, Standard Rectangular	2	Far UV Quartz	3.5	\$ 233.00
3-Q-10-MS	Stirring Cell, Fluorometer Rectangular	4	Far UV Quartz	3.5	333.00
9-Q-10-MS	Stirring Cell, Semi-Micro Rectangular	2	Far UV Quartz	1.8	307.00
18-Q-10-MS	Stirring Cell, Micro Rectangular	2	Far UV Quartz	0.9	333.00
21-Q-10-MS	Stirring Cell, Standard with Stopper	2	Far UV Quartz	3.5	299.00
23-Q-10-MS	Stirring Cell, Fluorometer with Stopper	4	Far UV Quartz	3.5	403.00
29-Q-10-MS	Stirring Cell, Semi-Micro with Stopper	2	Far UV Quartz	1.8	346.00
28-Q-10-MS	Stirring Cell, Micro with Stopper	2	Far UV Quartz	0.9	364.00
MSB-5X2	Magnetic Stirring Bars for above cells, 5mn	n x 2mm, 10/Pack			\$ 55.00/pack
MSB-6X1.5	Magnetic Stirring Bars for above cells, 6mn		\$ 55.00/pack		
MSB-6X3	Magnetic Stirring Bars for above cells, 6mn	n x 3mm, 10/Pack			\$ 55.00/pack

Magnetic Stirrer

Spinette Magnetic Cell Stirrer


The Spinette cell stirrer is a simple solution to mixing solutions inside of a spectrophotometer or fluorometer cell. The spinette drives a small electronic magnetic stirrer which has the same exterior dimension as a normal cell and is only 5mm thick. It is placed into the bottom of the cell compartment, the cell is then placed on top of the magnet drive. A ribbon connector is brought up the non-transmitting side of the cell and then connected into the spinette controller. A standard cell is fully mixed within 4 seconds. The spinette works well with any standard cell with a 10mm x 10mm interior or with any of the special stirring cells as listed above. Each spinette comes complete with one stir bar, electronic magnet driver, ribbon cable and 120V AC controller.

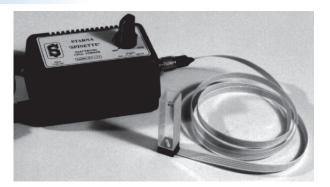


Catalog Number	Description	Price
9400	Spinette Cell Stirrer	\$ 979.00 ea.
MSB-5X2 MSB-6X1.5	Magnetic Stirring Bars, 5 x 2mm, 10/Pack Magnetic Stirring Bars, 6 x 1.5mm, 10/Pack	55.00/pk 55.00/pk
MSB-6X3	Magnetic Stirring Bars, 6 x 3mm, 10/Pack	55.00/pk
	www.starnacells.com	37

Constant Temperature Cells


Constant temperature cell Types 62, 63 and 65 are designed with a separate chamber surrounding the sample compartment through which a temperature controlling fluid can be circulated to keep the sample at a specific temperature. Type 62 cells have connection tubes exiting from the top, Type 63 on the front face and Type 65 on the upper portion of the cell. The light beam of the spectrophotometer is not intended to pass through the water jacket. PTFE stoppers are supplied as pictured.

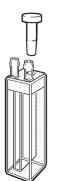

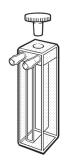

Type 63

Type 62

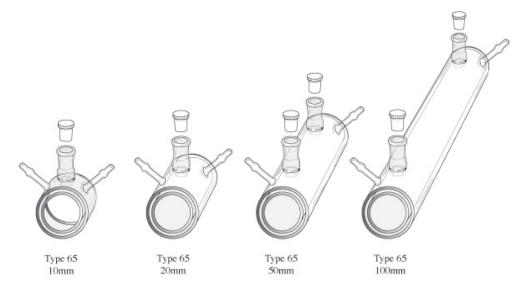

Catalog	Path Length	E	xterior, r	nm	Sample	e Cham	ber, mm	Overall	Nominal	Polished	Price
Number	mm	Width	Length	Height	Width	Lengtl	1 Height	Height	Vol ml	Windows	per cell
Spectrosil [®] F	ar UV Quartz,	windows	Useabl	range 1	70 to 270	0 nm			Rectangula	r, Vertical Fl	ow Tubes
62-Q-10	10	12.5	12.5	48	4	10	40	60	1.520	2	\$ 556.00
								Re	ectangular, H	Horizontal Fl	ow Tubes
63-Q-10	10	12.5	12.5	48	4	10	40	60	1.520	2	\$ 556.00

Catalog	Path Length	Exterior	, mm	Sample Cha	mber, mm	Overall	Nominal	Polished	Price
Number	mm	Diameter	Length	Diameter	Length	Height	Vol ml	Windows	per cell
Spectrosil [®] I	Far UV Quartz w	vindows, Usea	ble range: 1	70 to 2700 nm				(Cylindrical
65-Q-0.01	0.01	22	20	10	0.01	32	0.737	2	\$ 1,182.00
65-Q-0.1	0.1	22	20	10	0.1	32	0.747	2	973.00
65-Q-1	1	22	20	10	1	32	0.849	2	848.00
65-Q-2	2	22	20	10	2	32	0.962	2	830.00
65-Q-5	5	22	20	10	5	32	1.300	2	778.00
65-Q-10	10	22	20	13	10	32	0.825	2	391.00
65-Q-20	20	22	20	13	20	32	1.650	2	438.00
65-Q-50	50	22	50	13	50	32	4.125	2	587.00
65-Q-100	100	22	100	13	100	32	8.250	2	712.00

38

Cell Cleaning Suggestions

Why Clean Cells?

Clean cells are the foundation of any spectrophotometric or fluorometric analysis. The residue from previous analysis will cause inaccuracies, low sensitivity and lack of precision. More important, it will waste your time! Also inspect the condition of the cells. If they are cracked, chipped or scratched it is important to replace the cells with new ones as your time is more valuable than the cost of new cells.

How to Clean Cells: It is important to determine the residual material in the cell that needs to be removed. The table below will give you some suggestions for cleaning:

Solvent	Material	Suggested Cleaning Methods
Aqueous	Protein, DNA, Biologics	Warm water with detergent, Dilute acid rinse, Copious water rinse
Aqueous	Salt solution	Warm water Acid rinse, copious water rinse
Aqueous	Basic solutions	Warm water with detergent, Dilute acid rinse, Copious water rinse
Organic	Oil based	Rinse with solvent, Warm water with detergent, Dilute acid rinse, Copious water
rinse		
Organic	Alcohol solutions	Rinse with solvent, Copious water rinse
Organic	Acidic solutions	Rinse with solvent, Copious water rinse
Organic	Basic solutions	Rinse with solvent, Dilute acid rinse, Copious water rinse
Fluorescence mea	asurements - Clean cells in N	itric Acid (5M) use a copious water rinse immediately before use.

General Considerations:- Keeping cell clean while in use is the most important element of a long, useful Cell life. During the day, never let your cells dry out. If you keep them in a water or solvent bath between usage, the material that you are using will not have a chance to dry out and stick. Use only lens cleaning paper or fine cloth to wipe the optical surfaces, most paper products contain wood fibres which may scratch or damage the cell face or surface. At the end of the day, ensure all cells are well cleaned and store in a suitable container after drying.

Definitions

Dilute Acid	Dilute Hydrochloric acid (2M) or Nitric Acid (2M)
Acid	Hydrochloric (5M) acid or Nitric Acid (5M) (see note below)
Solvent rinse	Rinse with the solvent that solvated your sample in the first place!
Copious water rinse	Use a pure water like deionized, distilled or RO and rinse at least 10 times
Detergent	Use a neutral pH detergent if available but dilute acid wash and water rinse to remove detergent residues

Important Exceptions:


5M Nitric acid: Do not use this treatment on Anti-reflection or mirror coated cells Ultrasonic Cleaners: We do not recommend the use of ultrasonic cleaning baths with cells. Each bath generates a different frequency and if your bath operates at the resonant frequency of a cell, the cell will break. We do not warranty our cells for cleaning in an ultrasonic cleaner.

Cell Holder for Infrared Instruments

Cell Holder, long path length

Cell holders Type CH-1/50 and CH-1/100 are for rectangular cells up to 50 and 100mm path lengths respectively. All come complete with a 3" x 2" back plate to fit standard Infrared instrument holders.

Cat. No	Description	Price, each
CH-1/50	Rectangular to 50mm	\$ 372.00
CH-1/100	Rectangular to 100mm	456.00

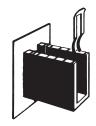

Starna Cells, Inc. PO Box 1919 Atascadero, CA 93423

Terms of Sale

Our terms of sale are Mastercard/VISA, American Express or net 30 days to recognized accounts. If you are unsure of your account status, please call us so that we can verify your account for you. The prices published in this catalog are current at the time of publication and are subject to change without notice.

Our terms of shipment are FOB Atascadero, CA USA.

Method of Shipment

We usually ship via United Parcel Service and have all levels of service available including Red (overnight). If required we can also ship via Federal Express. Foreign orders are usually shipped via Federal Express or UPS. All shipping charges are prepaid and added to the invoice.

In Stock

A great effort is made to keep most items in stock so that your order may be shipped on the same day that we receive it. However, at times we are out of stock and we will backorder any products that we could not ship to you immediately. Your backordered items will ship on the same day that they are received from our manufacturing facility.

Cell Matching

With common cells, such as the standard rectangular cell type 1-Q-10, the size of the matched set is virtually unlimited. The less commonly used cells can usually be matched in sets of 2 or 4 cells. All cells are priced as single cells but, when you order multiples of the same cell, we make every effort to match the cells as closely as possible.

Product Warranty

Our Starna[®] brand Spectrophotometer and Fluorometer cells are warranted to meet the specifications for transmission as listed in the specification table on page two of this catalog. In addition, the dimensions of each cell are warranted to be within the description of physical dimensions listed in this catalog. If the transmission or dimensions are not to your satisfaction, please call us within fifteen days of your receipt of the cells for warranty repair or replacement. We perform extensive quality control over our production and use only branded materials to assure you receive a cell that will perform to the maximum specification of your instrument. Any goods to be returned under warranty require a Return of Merchandise Authorization (RMA) number. The RMA can be obtained by calling our customer service department.

Starna Cells, Inc

Mail orders:

Starna Cells, Inc PO Box 1919 Atascadero, CA 93423

Phone Orders: (800) 228-4482 or

(805) 466-8855

E-Mail Address:

sales@starnacells.com

FAX Orders: (805) 461-1575

World Wide Web Site:

http://www.starnacells.com

Technical Information

We have technical staff available to assist you with any information that you may need for the selection of a cell. If you are unsure of the material or physical configuration of a cell, please give us a call and we will be pleased to answer any questions that you may have. For advice on the cleaning of cells please see the cleaning guide on page 39 of this catalog.

Our hours of business are 8:30AM to 4:30PM, Pacific Time